


pNet

(Personalized Functional Network Modeling)

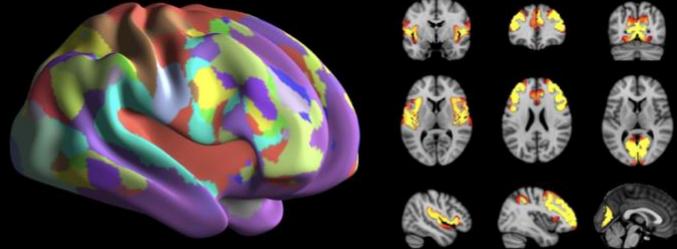

Help Document


A guideline for pNet

Yuncong Ma

Toolbox Overview

This toolbox is designed to provide a user-friendly interface to perform personalized functional network computation. It is open-source, cross-platform, and expandable, making it easy to use for users with little knowledge about computation, customizable for users with sufficient coding skills. It allows flexible data search for fMRI scans with support for different data formats. To initialize the computation of personalized functional networks (FN), this toolbox allows user to compute a new group-level FN based on the data available, or load a precomputed group-level FN using this toolbox, or from other methods, including group-ICA and GIGICA. Flexible visualization options enhance fast result navigation and customizable in-depth comparisons. It also provides several simple statistical analysis methods to investigate the correlation between the topology of personalized FNs to behavior data. In addition, quality control is integrated to ensure the spatial correspondence of personalized FNs with respect to the corresponding group-level FNs.

This toolbox can be downloaded from <https://github.com/YuncongMa/pNet>

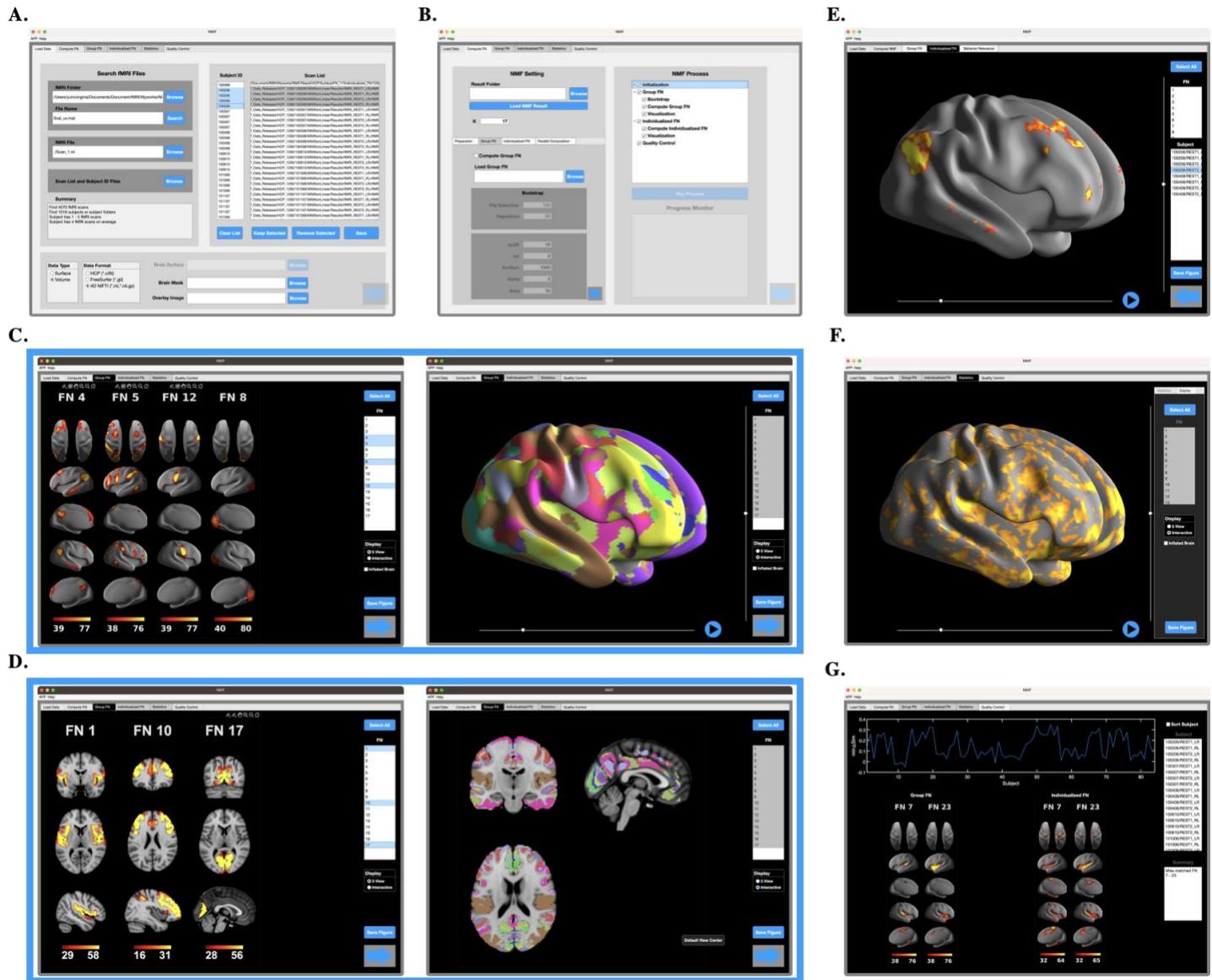


pNet app snapshots

- (A) Welcome page for the toolbox.
- (B) A module for loading fMRI scans and brain template files.
- (C) A module to setup computation parameters for both the group-level and individualized FNs.
- (D) Surface-based visualization of both group and personalized FNs ($k=17$) using HCP S1200 dataset, with left panel showing a binarized atlas generated from the group FNs and the right panel showing five views of three personalized FNs. All the color bars of intensity maps were set from the maximum value of the map to its half value.
- (E) Volume-based visualization of both group and personalized FNs ($k=17$) from Zhen's multi-cohort iSTAGING study, with left panel showing a binarized functional atlas generated from the group FNs and the right panel showing a three-slice view of three personalized FNs.
- (F) Surface-based visualization of the maximum t value (two sample t-test, FDR correction, p-


value=0.001) of sex differences of individualized FNs (k=50) of the HCP S1200 dataset.

(G) A module for quality control, showing one scan with two FNs mismatched to their group-level counterparts.

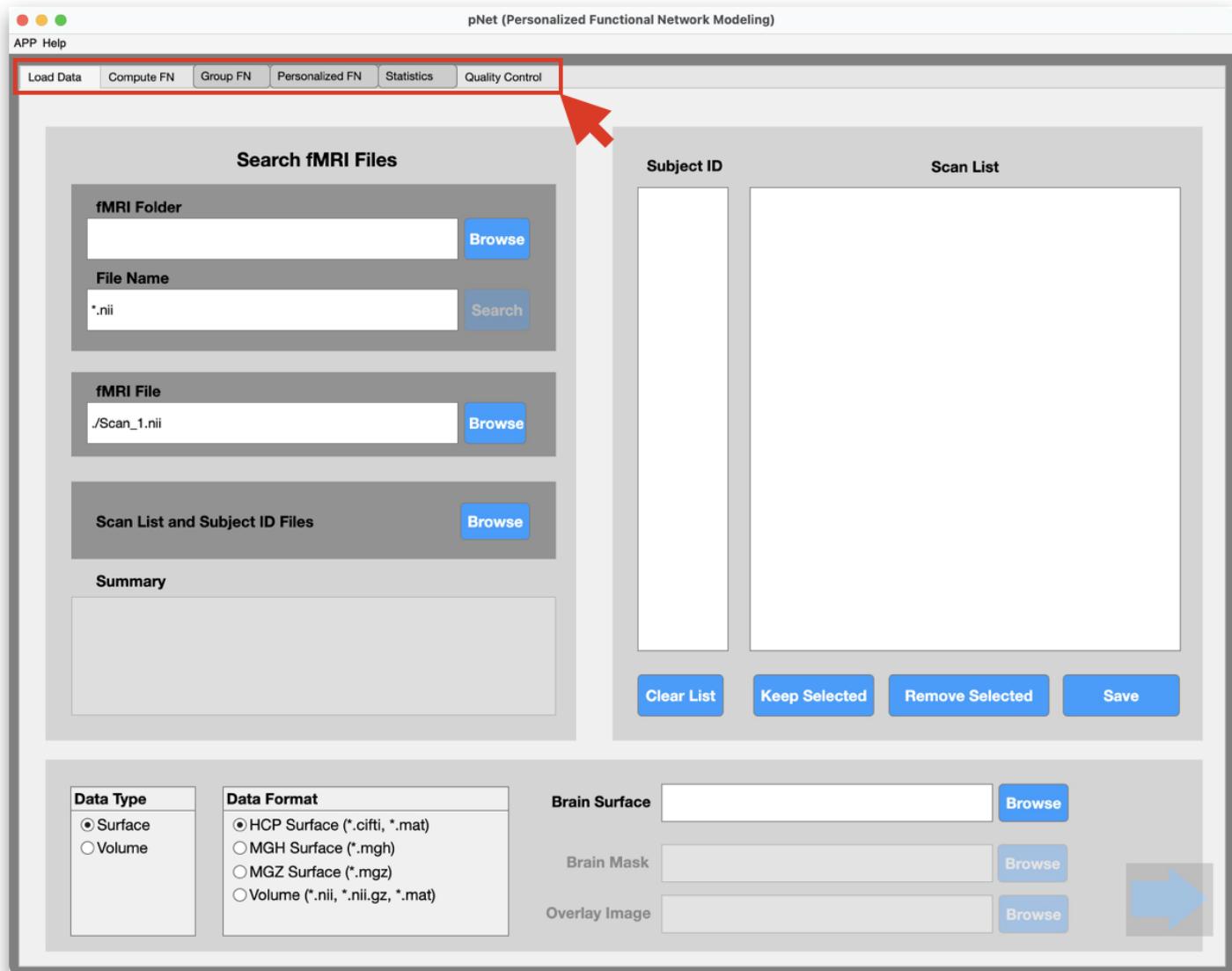


Tabs (highlighted in the red box) for different modules to search data, setup FN computation, check results of group-level and personalized FN, conduct statistics, and quality control.

pNet provides a unified and intuitive design of UI and data structure. It consists of six modules

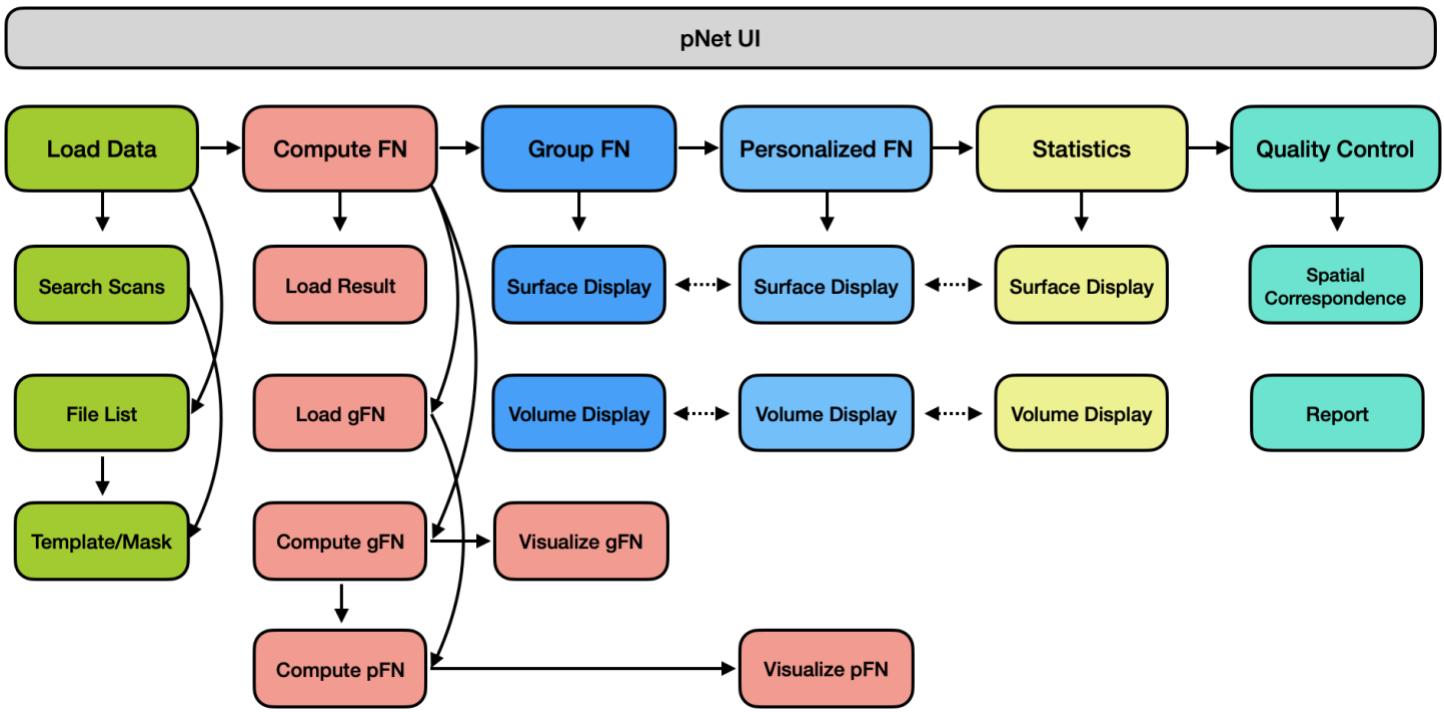

The workflow of pNet UI.

It consists of six major modules, including data loading, FN computation, group-level FN display, personalized FN display, statistics, and quality control. The six modules correspond to six major UI panels in pNet, and six subfolders in the pNet result folder. The pNet UI is designed to follow the left-to-right workflow direction to setup each step. Curved arrows show alternative workflows. The FN computation module (colored in light red) can use precomputed group-level FNs to derive personalized FNs as one option. Precomputed visualization of group-level and personalized FNs will be used in the two display modules (group FN and personalized FN). pNet features synchronized display to facilitate the comparisons between group-level and personalized FNs, as well as statistical results. The dashed double-sided arrows denote this display synchronization between three modules.

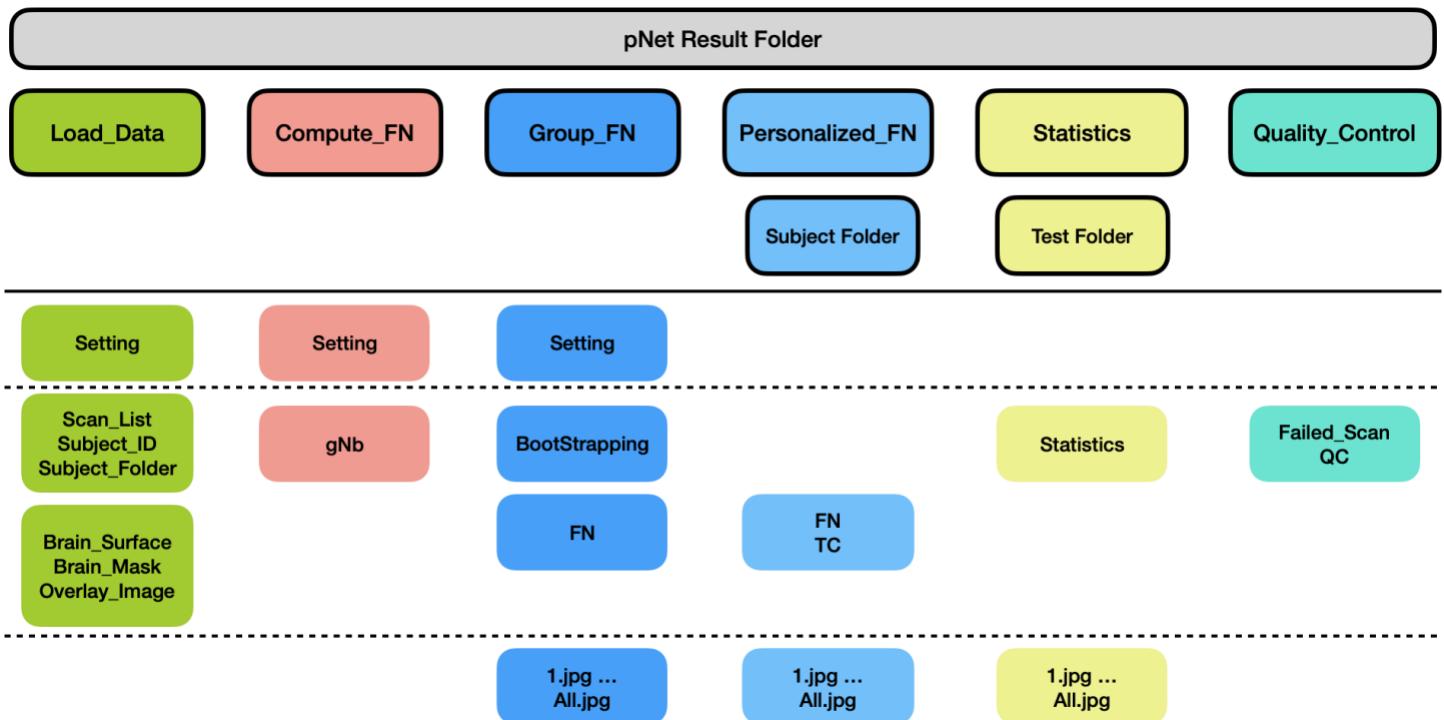

Data structure in the pNet result folder.

Folders are noted in boxes with black borders above the black horizontal line, and files are in boxes without black border, separated into three categories (setting, data, and figure files) by the dashed lines. Boxes are color coded for different modules. In the first three modules, the setting files are present. In the personalized FN folder, subfolders are named by subject information. In statistics folder, subfolders are named by their method information.

Content

APP INSTALLATION	9
SYSTEM AND HARDWARE REQUIREMENT.....	9
RUN IN MATLAB	10
INSTALL PNET AS A MATLAB APP.....	10
RUN COMPILED PNET WITH MATLAB RUNTIME.....	13
COMPILE PNET AS A STANDALONE SOFTWARE	13
DATA LOADING.....	15
SEARCH FOR fMRI FILES.....	15
1. <i>Search fMRI files based on their extension in a folder.</i>	16
2. <i>Select fMRI files</i>	19
3. <i>Load prepared file list information in txt format.</i>	20
ORGANIZED fMRI SCAN FILES	21
1. <i>Clear, select or remove fMRI scan files</i>	21
2. <i>Save the fMRI scan list and subject ID information into text files</i>	22
LOAD BRAIN TEMPLATE.....	22
1. <i>Surface data type</i>	22
2. <i>Volume data type.</i>	23
SETUP FN COMPUTATION	25
SETUP RESULT FOLDER.....	25
LOAD PRE-COMPUTED RESULTS	25
SETUP GROUP-LEVEL FN	26
1. <i>Setup the group-level FN computation.</i>	27
2. <i>Load a precomputed group-level FN file</i>	28
SETUP PERSONALIZED FN COMPUTATION.....	29
SETUP PARALLEL COMPUTATION	30
RUN THE FN COMPUTATION.....	30
CHECK GROUP-LEVEL FN	32
SURFACE DATA	32
1. <i>Precomputed visualization</i>	32
2. <i>Interactive visualization</i>	34
VOLUME DATA.....	36
1. <i>Precomputed visualization</i>	36
2. <i>Interactive visualization</i>	38
CHECK PERSONALIZED FN	40

SURFACE DATA	40
1. <i>Precomputed visualization</i>	40
2. <i>Interactive visualization</i>	42
VOLUME DATA.....	44
1. <i>Precomputed visualization</i>	44
2. <i>Interactive visualization</i>	46
COMPARE GROUP-LEVEL AND PERSONALIZED FN.....	48
STATISTICAL ANALYSIS	49
ONE-SAMPLE T-TEST.....	49
TWO-SAMPLE T-TEST	52
QUALITY CONTROL.....	54
EXAMPLE	56

App Installation

This chapter provides guidelines about ways to run this toolbox. pNet only runs successfully on MATLAB no older than 2021a.

Users can use the toolbox in several ways as shown below:

1. Run from the MATLAB command window.
2. Install pNet as a MATLAB app.

System and hardware requirement

It is recommended to have at least 4 CPU cores, 16GB memory (RAM), and 10GB free disk space.


Please be aware that more CPU, memory and disk space will be needed to process bigger data or use parallel computation. MATLAB older than 2021a does not have some essential APP UI functions to give correct UI layout. As pNet is mostly tested on 2022b, it is recommended to use this version or later.

MATLAB	macOS (Intel CPU)	macOS (Apple Silicon)	Linux	Windows 10
2020b				✗
2021a	✓	✓	✓	✓
2021b	✓			✓
2022a	✓			✓
2022b	✓			✓
2023a	✓			✓



Run in MATLAB

This method requires to download all the functions and app icons to run pNet in MATLAB. Make sure that the subfolder Private_Function and the help document are available in the app folder.

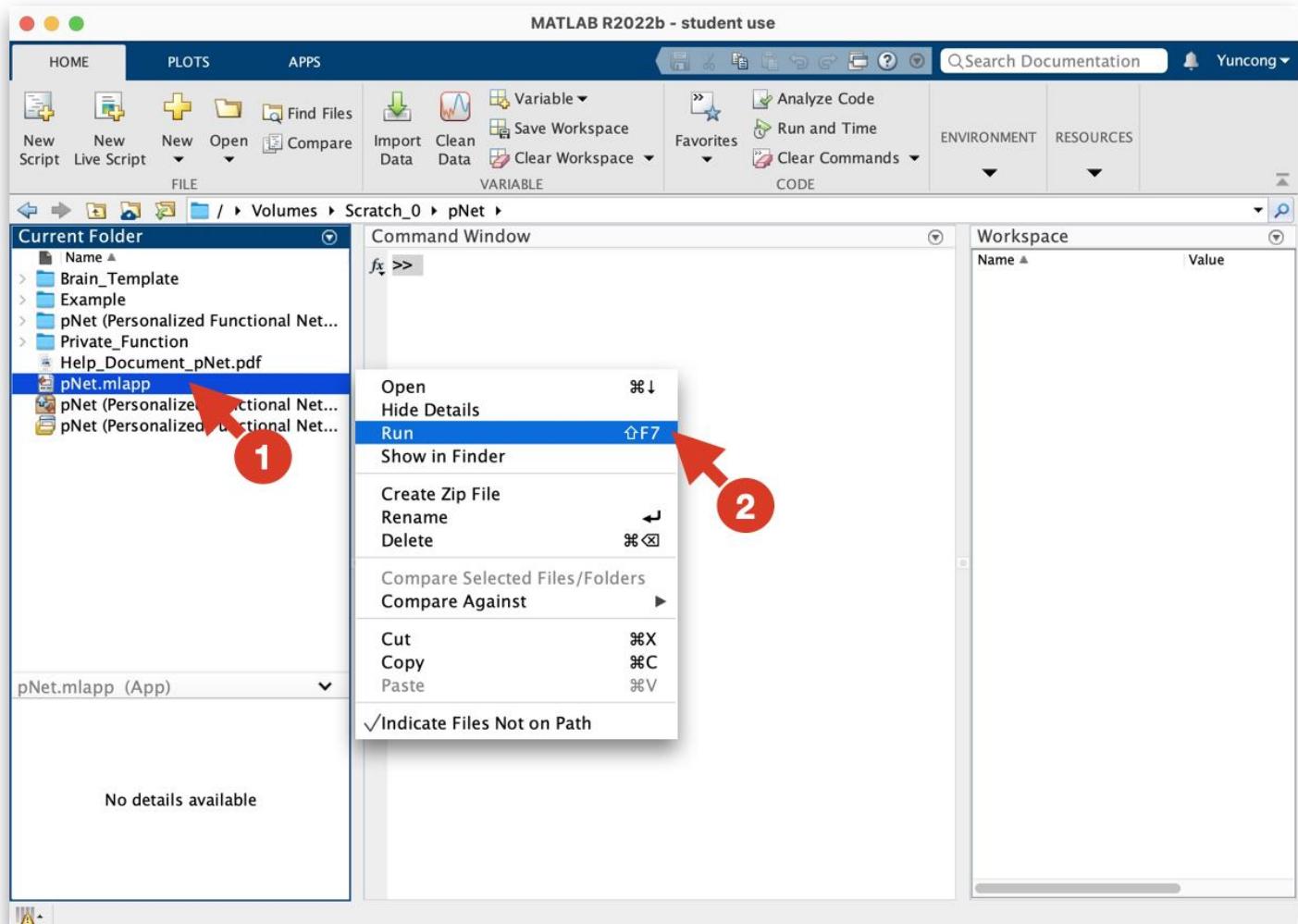


Run the pNet app.

(1-2) Open the pNet folder, and right click the pNet.mlapp, and click the run option.


Install pNet as a MATLAB app.

This way will allow users to use the app without navigating to the app folder, or adding the app folder to the path.

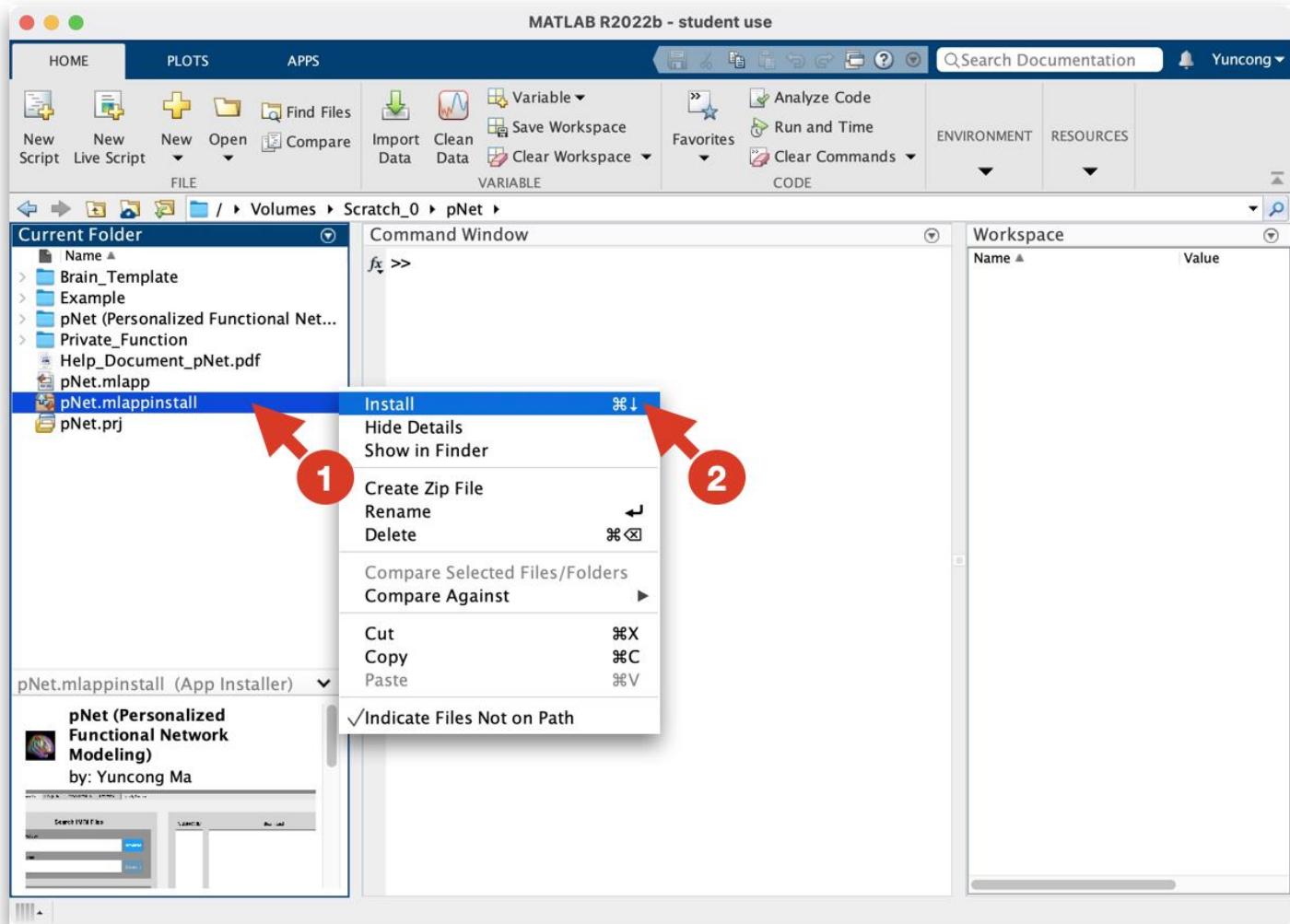

Install the pNet app UI.


(1-2) Double click or right click on the file “pNet.mlappinstall” to install the software on MATLAB.

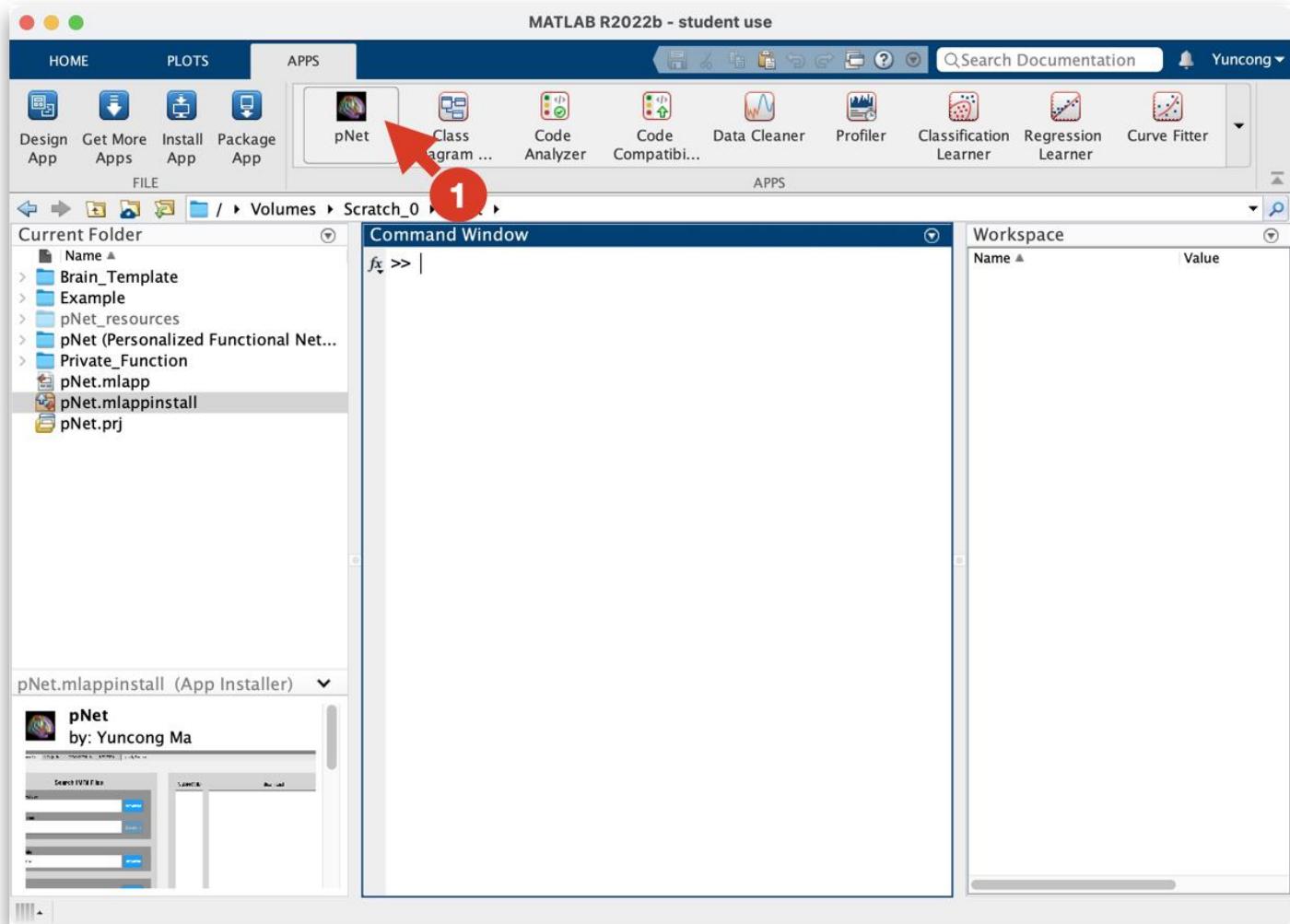

Run the pNet app from the APPS menu.

- (1) Open the APPS menu tab, and click on the pNet app.

Users will see a welcome page, and a main user interface (UI). The main UI contains multiple tabs on the top for different modules to setup the computation and visualization. Users could simply use the order from left to right to follow the workflow of this toolbox.

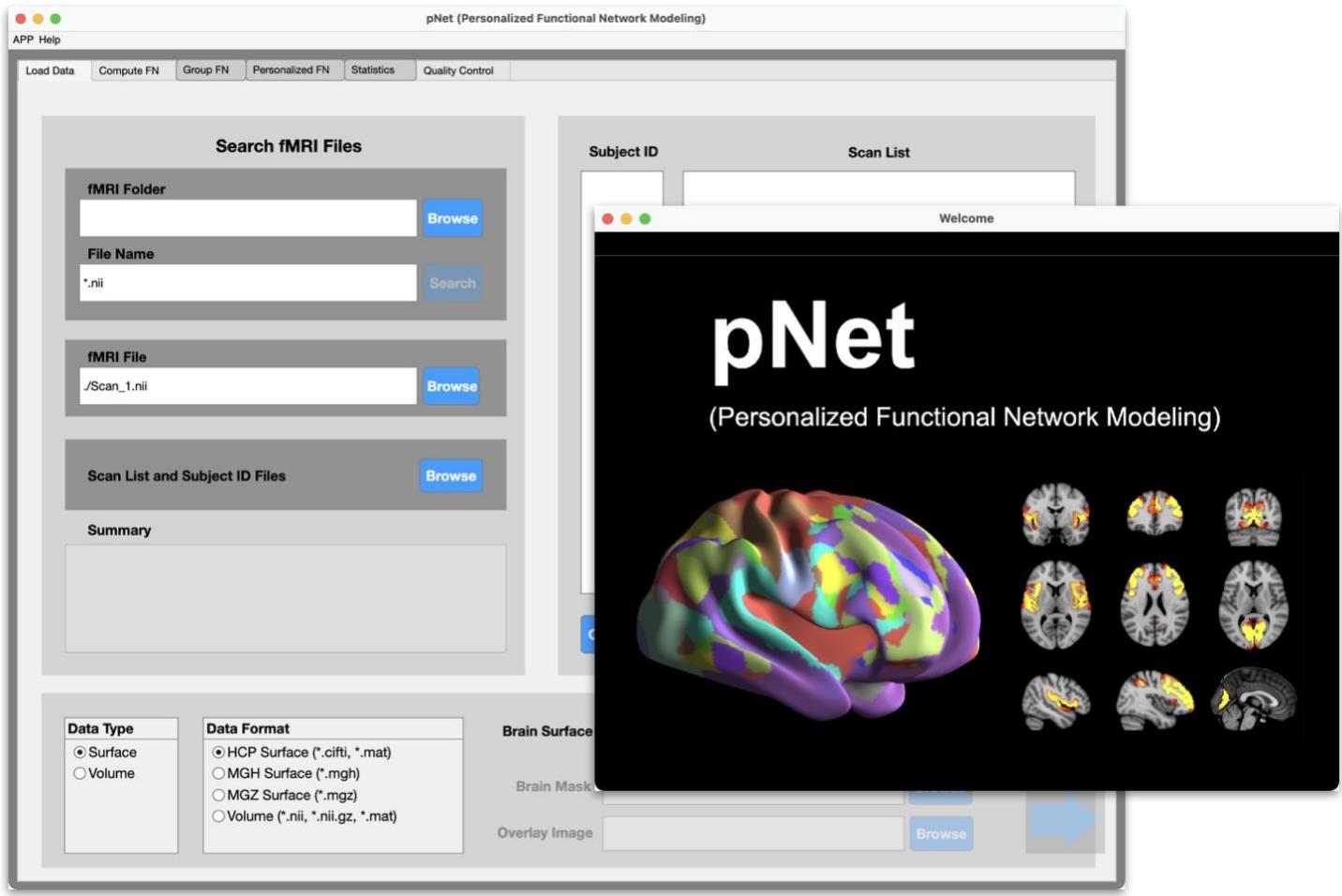

The app welcome page (right window with black background) and the main UI.

Run compiled pNet with MATLAB Runtime

There are two options to get a standalone software with or without the free MATLAB runtime software (<https://www.mathworks.com/products/compiler/matlab-runtime.html>). Several pre-compiled pNet versions are provided in the subfolder “Compiled_Runtime”. Download the version compatible to your computer’s operation system, and download and install the MATLAB runtime software. Double click the compiled pNet file (pNet.exe in Windows) to run the app. Details can be found in the README.txt in each subfolder containing the compiled pNet software.

Compile pNet as a standalone software

MATLAB (academic or commercial license) allows users to compile the source code with their modifications or preference of different MATLAB Runtime versions. Details can be found on MATLAB website: <https://www.mathworks.com/help/compiler/create-standalone-app-using-application-compiler-app.html>. When a MATLAB APP is deployed, some functions cannot be used, such as “addpath”. Also, the code will be operation-system dependent, if terminal command or other non-

MATLAB functions were used. This requires sufficient software development experience to ensure the compatibility. Some examples can also be found in the source code of pNet. Please check the startup function in the file pNet.mlapp.

Data Loading

This chapter provides guidelines to search and organize fMRI can files, load brain template files.

Search for fMRI files

Users need to firstly select the data type and data format.

The screenshot shows the 'Search fMRI Files' window of the pNet software. The window is divided into several sections:


- Search fMRI Files:** This section contains three input fields:
 - fMRI Folder:** A text input field with a 'Browse' button.
 - File Name:** A text input field containing '.nii' with a 'Search' button.
 - fMRI File:** A text input field containing '/Scan_1.nii' with a 'Browse' button.
- Scan List and Subject ID Files:** A section with a 'Browse' button.
- Summary:** A large empty text area.
- Subject ID:** A large empty text area.
- Scan List:** A large empty text area.
- Buttons:** At the bottom are four buttons: 'Clear List', 'Keep Selected', 'Remove Selected', and 'Save'.
- Bottom Section:** This section includes:
 - Data Type:** Radio buttons for 'Surface' (selected) and 'Volume'.
 - Data Format:** Radio buttons for 'HCP Surface (*.cifti, *.mat)' (selected), 'MGH Surface (*.mgh)', 'MGZ Surface (*.mgz)', and 'Volume (*.nii, *.nii.gz, *.mat)'.
 - Brain Surface:** A text input field with a 'Browse' button.
 - Brain Mask:** A text input field with a 'Browse' button.
 - Overlay Image:** A text input field with a 'Browse' button.
 - A large blue arrow pointing right.

Select the data type and format.

- (1) Select the data type.
- (2) Select the data format.

The toolbox offers three ways to find fMRI scan files.

1. Search fMRI files based on their extension in a folder.

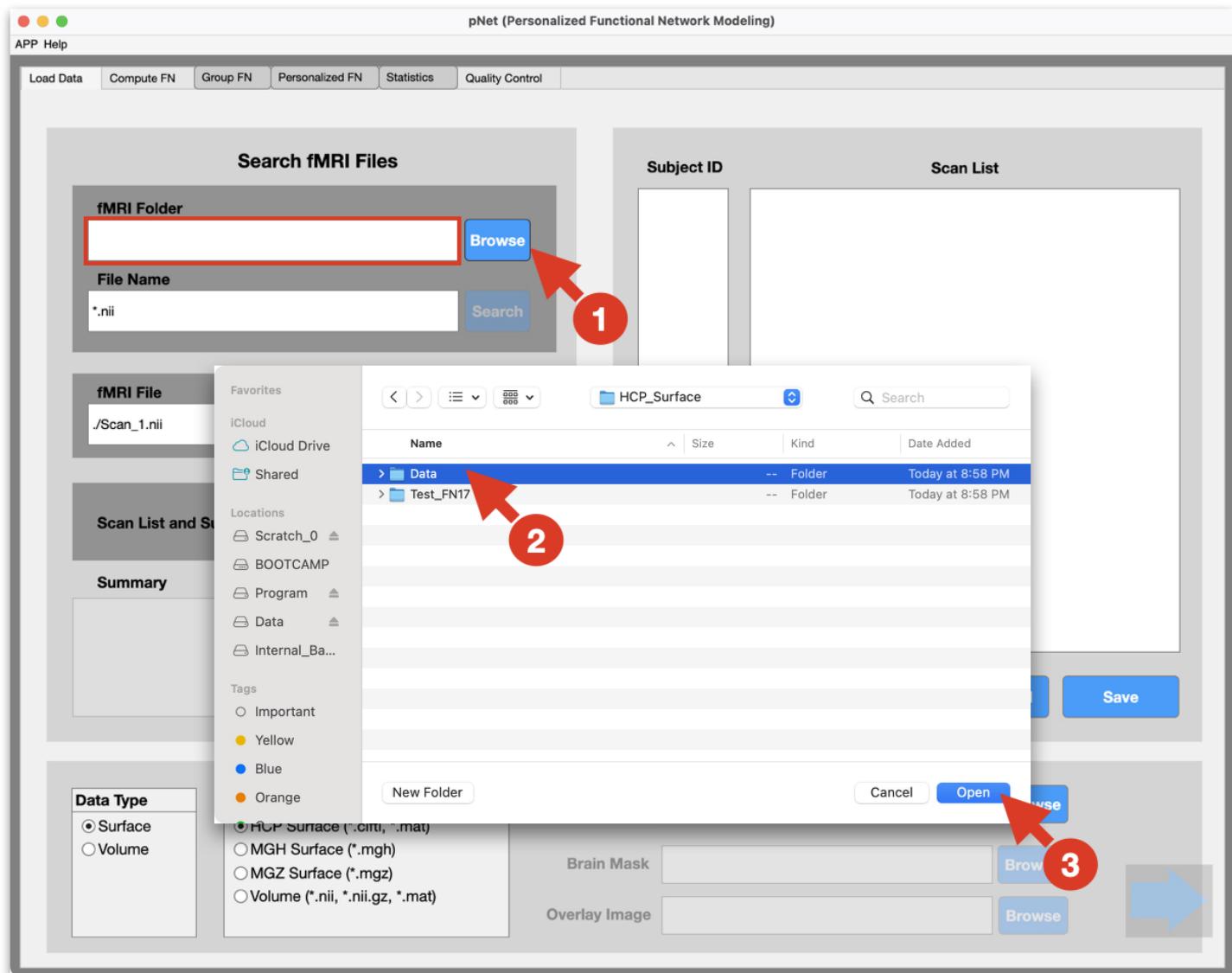


Select the data folder.

(1) Click browse button to open the folder selection window.

(2) Select the folder, and click the open button

(3) to finish the folder selection. The full folder directory information will appear in the white input field (highlighted in the red box). Users can also type in the full folder directory into the white input field (highlighted in the red box).

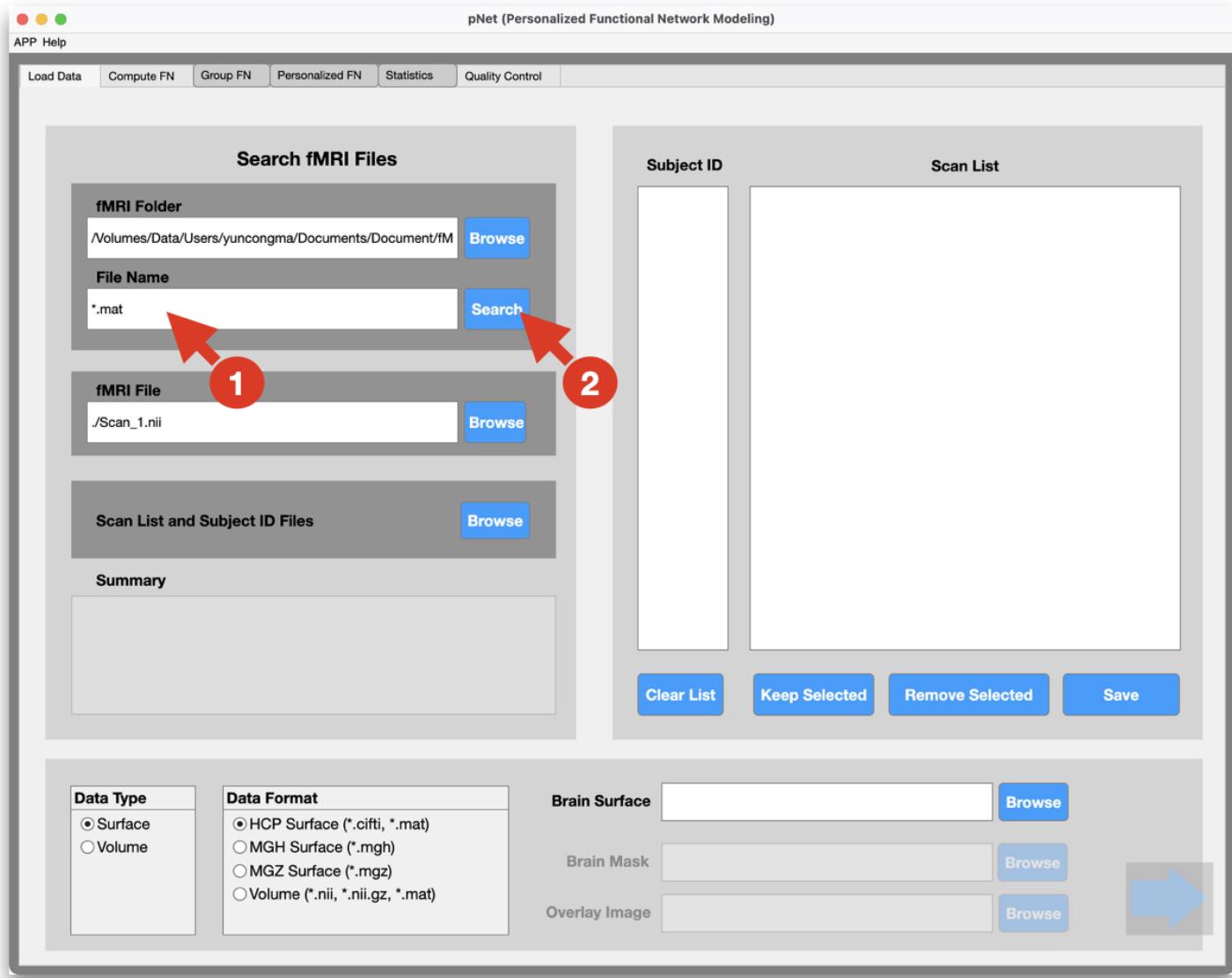

Setup the file extension information.

- (1) Specify the file extension, such as `*.mat` or `*.nii.gz`.
- (2) Click search button to perform file searching function.

pNet (Personalized Functional Network Modeling)

APP Help

Load Data Compute FN Group FN Personalized FN Statistics Quality Control

Search fMRI Files

fMRI Folder
 Browse

File Name
 Search

fMRI File
 Browse

Scan List and Subject ID Files **Browse**

Summary

Find 20 fMRI scans
 Find 10 subjects or subject folders
 Subject has 2 - 2 fMRI scans
 Subject has 2 fMRI scans on average

Subject ID	Scan List
100206	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
100206	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
100307	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
100307	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
100408	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
100408	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
100610	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
100610	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
101006	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
101006	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
101107	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
101107	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
101309	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
101309	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
101410	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
101410	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
101915	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
101915	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
102008	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res
102008	/Volumes/Data/Users/uncongma/Documents/Document/fMRI/Myworks/NMF/Res

Clear List **Keep Selected** **Remove Selected** **Save**


Data Type
 Surface
 Volume

Data Format
 HCP Surface (*.cifti, *.mat)
 MGH Surface (*.mgh)
 MGZ Surface (*.mgz)
 Volume (*.nii, *.nii.gz, *.mat)

Brain Surface **Browse**
Brain Mask **Browse**
Overlay Image **Browse**

Check the summary and navigate through the extracted subject ID and corresponding file directory. Summary module highlighted in the red box shows the total number of files, the range of file number per subject, average file number per subject. Error messages will also be shown if there is any.


2. Select fMRI files

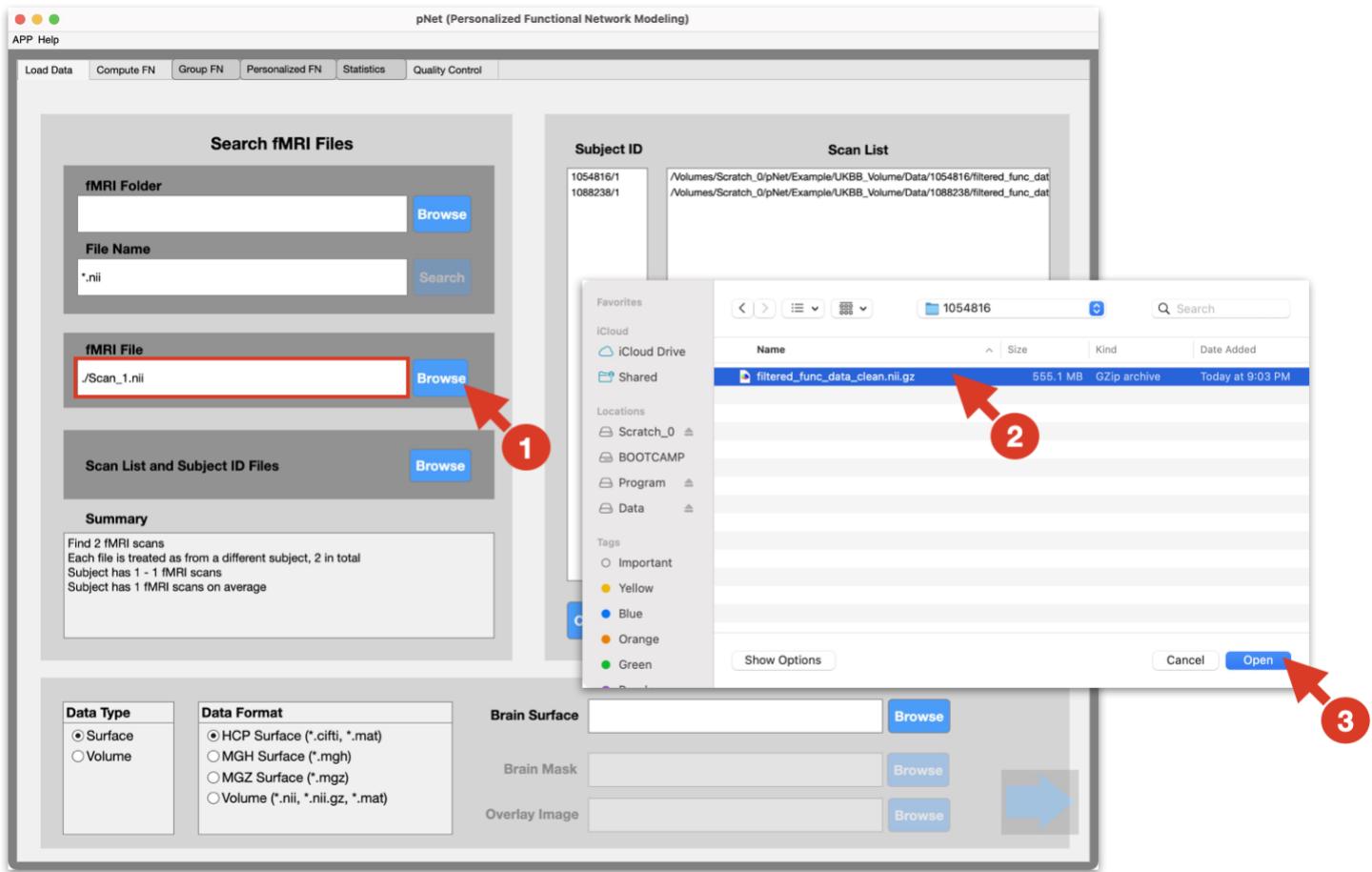

Select fMRI files.

- (1) Click the browse button.
- (2) Manually select one or multiple fMRI files with supported extensions, and click the open button
- (3) to load the full file directories. Users can also use the input field (highlighted in the red box) to type in the full file directories.

3. Load prepared file list information in txt format.

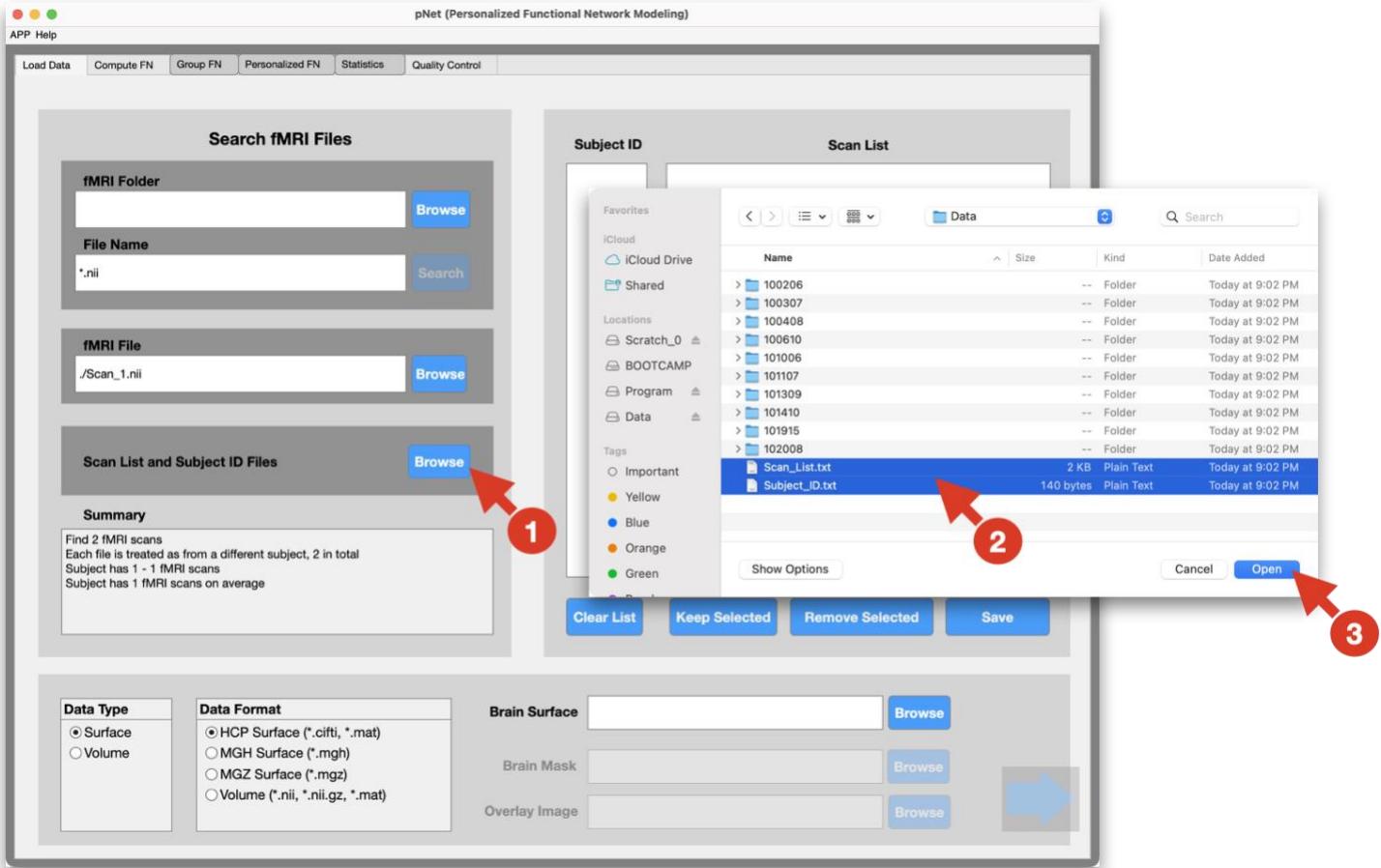


Load fMRI scan list and subject ID files.

- (1) Click the browse button.
- (2) Choose the text file containing full directories of fMRI scan files, and another text file containing subject ID information for each corresponding fMRI scan. The subject ID information file is not necessary. When it is not provided, automatic subject ID information will be extracted.
- (3) Click the open button to read those files.

Organized fMRI scan files


1. Clear, select or remove fMRI scan files

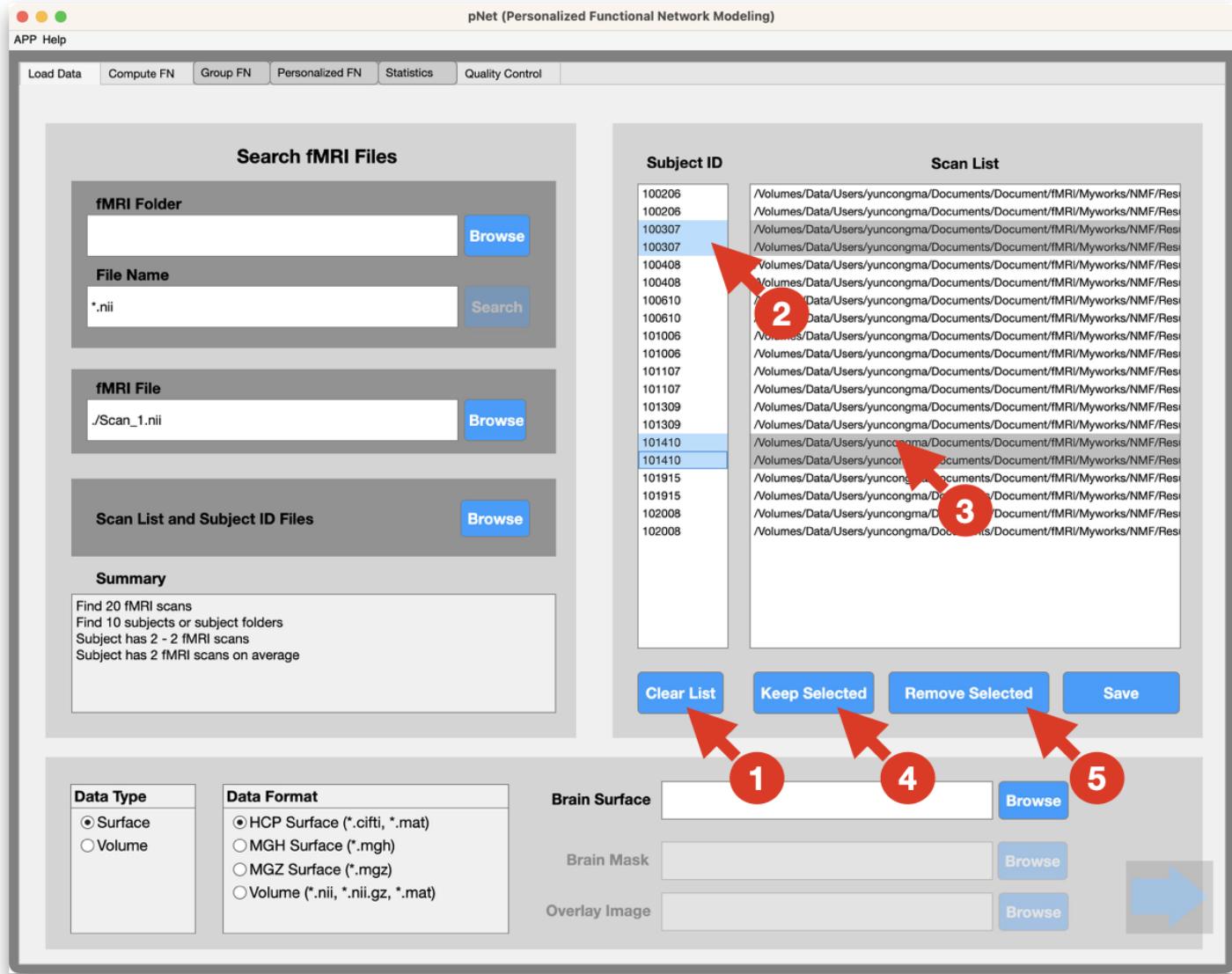

Remove, select or remove fMRI scan files.

- (1) Click the clear list button to empty the subject ID and scan list.
- (2-3) Select one or multiple files (hold command key in macOS, or hold control key in Windows and Linux) based on their subject ID or scan list.
- (4) Keep the selected fMRI scan files.
- (5) Remove selected fMRI scan files.

2. Save the fMRI scan list and subject ID information into text files

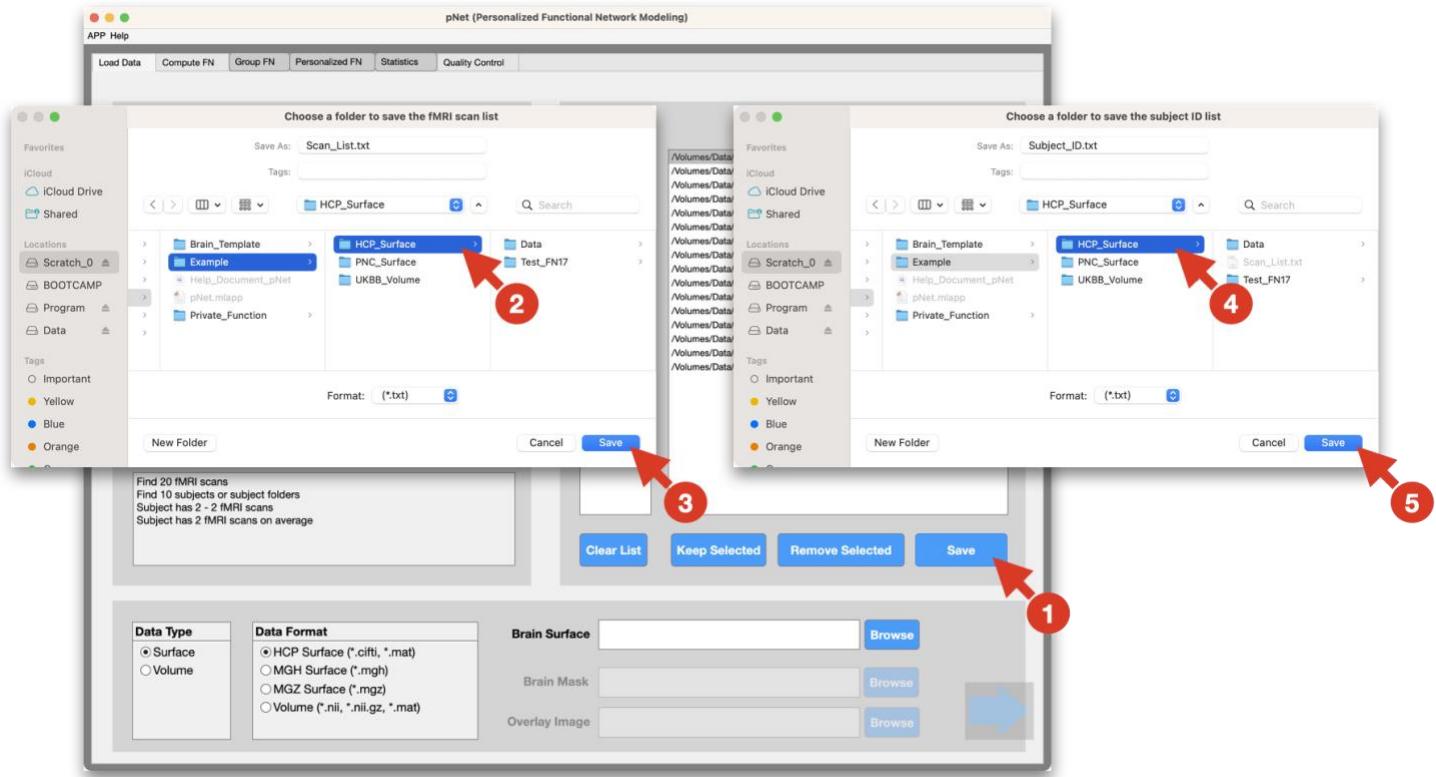

Save the subject ID and scan list information into text files.


- (1) Click the save button. A new window will be promoted for saving the scan list.
- (2) Select the folder to store the scan list file with the name “Scan_List.txt” in default.
- (3) Click the save button to save the scan list file.
- (4) Select the folder to store the scan list file with the name “Subject_ID.txt” in default.
- (5) Click the save button to save the subject ID file.

Load brain template

Surface and volume data types require different brain template files for both computation and visualization of group-level and personalized FNs.

1. Surface data type.

It requires to load a file called “Brain_Surface.mat”. This file a variable named as “Brain_Surface” with three fields: “Shape”, “Shape_Inflated”, “MW”.

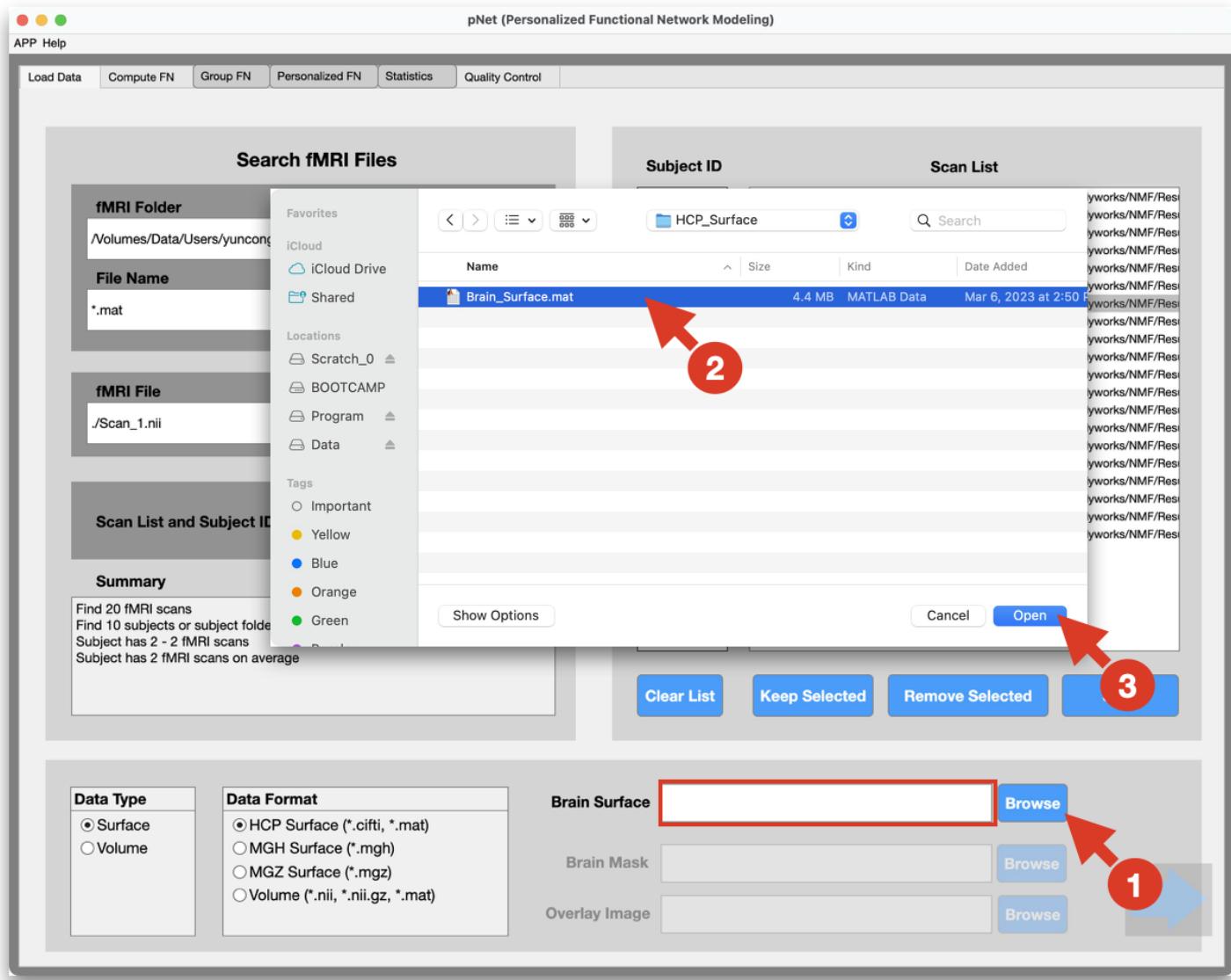


Load the brain surface file for surface data type.

- (1) Click the browse button.
- (2) Select the file “Brain_Surface.mat”.
- (3) Click open or double click the file to load it. Full file directory of the brain surface file will appear in the input field (highlighted in the red box). Users can also type in this information directly into the input field.

2. Volume data type.

It requires to load two files: “Brain_Mask.mat” and “Overlay_Image.mat”. The brain mask file is a binarized 3D mask to specify voxels to contain useful fMRI signals. The overlay image is a 3D image file that will be used as a background for overlaying group-level and personalized FNs.

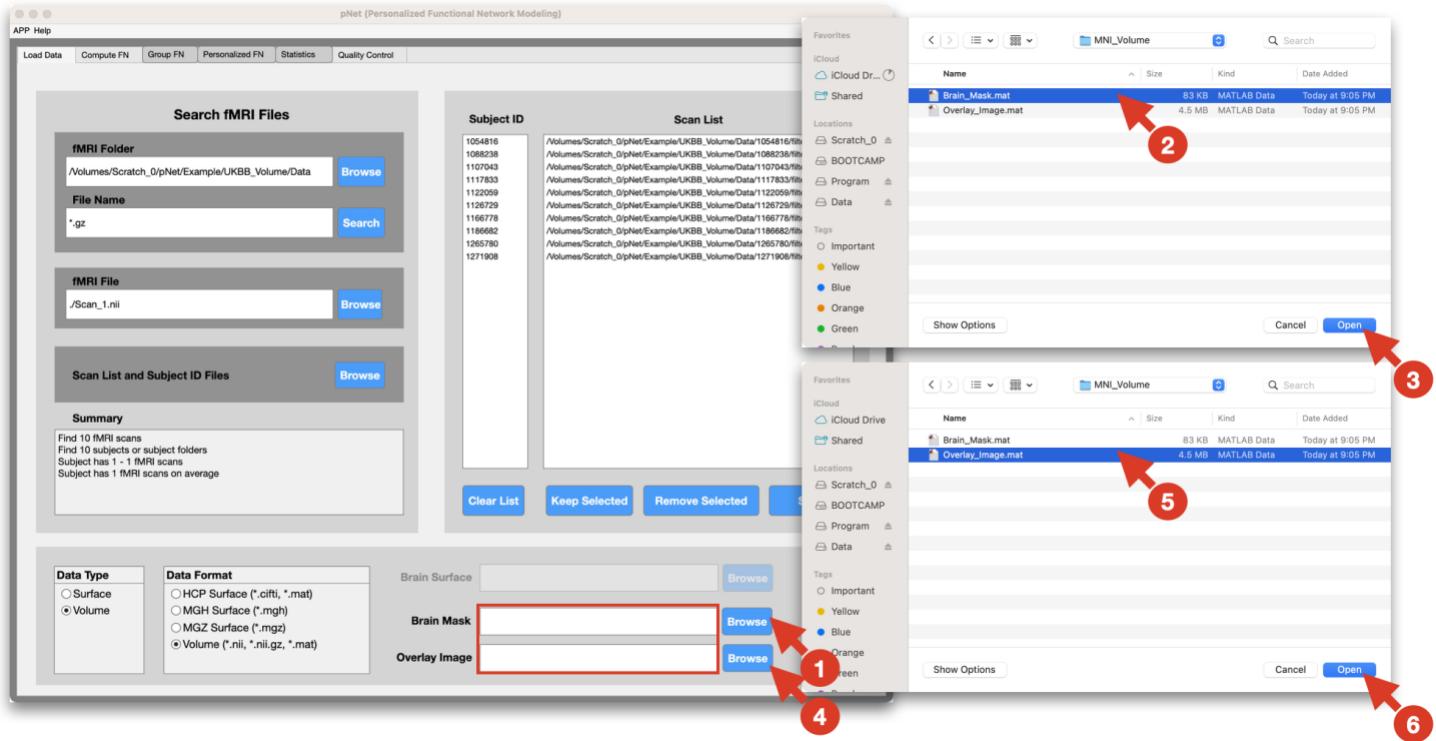


Load brain templates for the volume data type.

- (1) Click the browse button to load brain mask file, or use the input field (highlighted in the red box) to type in the full directory of the brain mask file. Example files can be found in the toolbox subfolder “/Brain_Template/MNI_Volume/Brain_Mask.mat”.
- (2-3) Select the brain mask file, and click open button to load.
- (4) Click the browse button to load the overlay image file, or use the input field (highlighted in the red box) to type in the full directory of the overlay image file. Example files can be found in the toolbox subfolder “/Brain_Template/MNI_Volume/Overlay_Image.mat”.
- (5-6) Select the overlay image file, and click open button to load.

Setup FN Computation

Setup result folder

The toolbox store all the fMRI file information, brain template files, results and figures of group-level and personalized FN, and other files in one result folder. The fMRI files will not be stored in the result folder to save storage space.

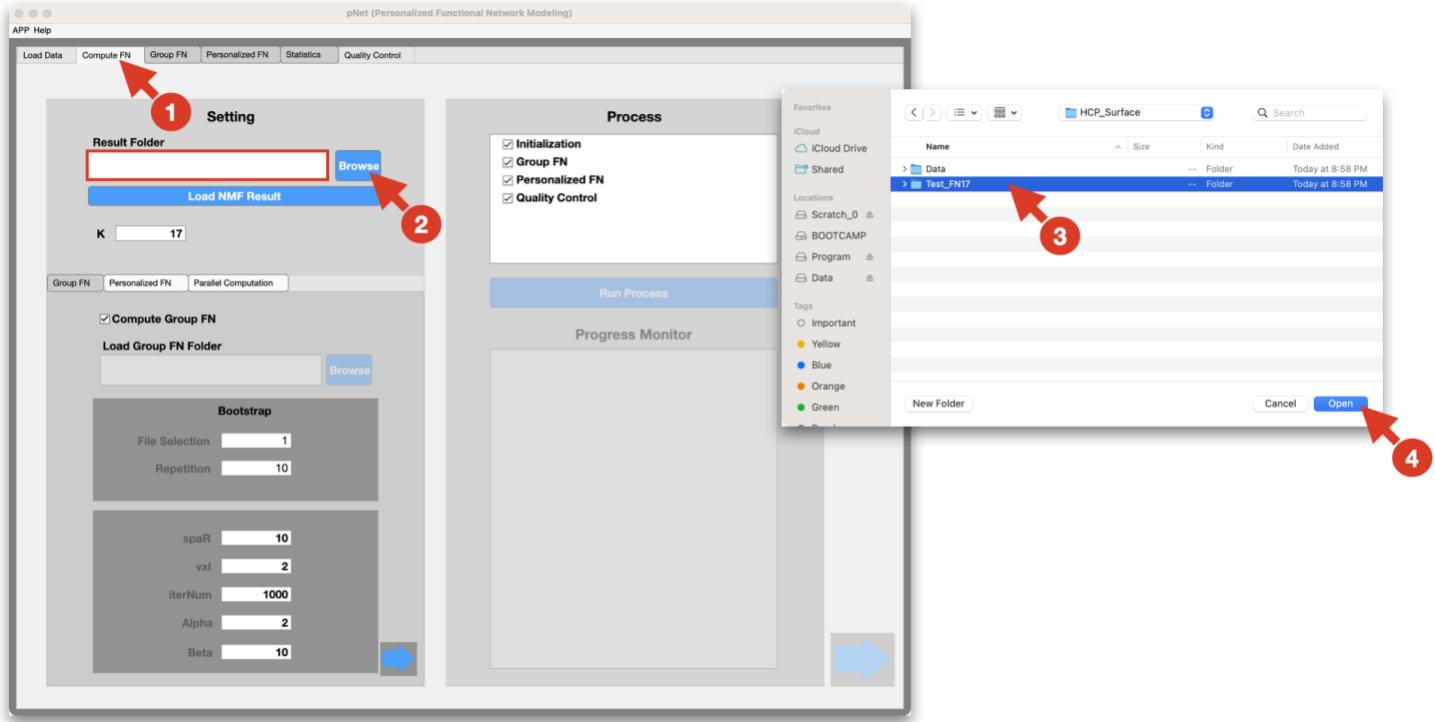


Select the result folder.

- (1) Navigate to the module “Compute FN”.
- (2) Click the browse button to select the result folder, or type in the full directory of the desired result folder.
- (3-4) Select the folder to store result, and click open button.

Load pre-computed results

The toolbox can automatically load pre-computed results for convenience. Users just need to simply select the result folder, all the information including fMRI scan file and subject ID, number of FNs, parameters of personalized FN computation will appear on the user interface.

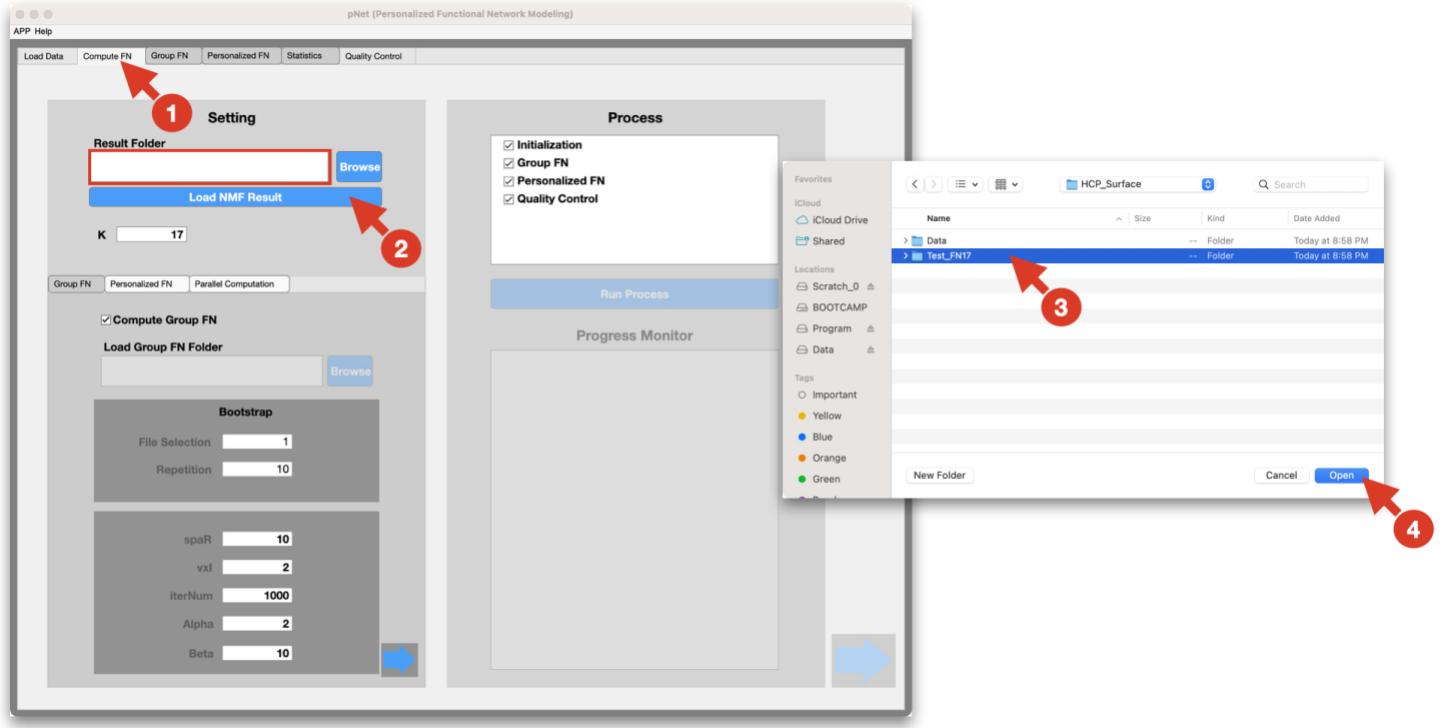


Select the folder containing precomputed results.

- (1) Navigate to the module “Compute FN”.
- (2) Click the browse button to select the result folder, or type in the full directory of the desired result folder (highlighted in the red box).
- (3-4) Select the folder to store result, and click open button. The toolbox will automatically check the integrity of the result folder. Then it will load results and jump to the group FN module to show precomputed group FN figures.

Setup group-level FN

The toolbox allows users to either load a precomputed group-level FN or perform the group-level FN computation using the data available.


1. Setup the group-level FN computation

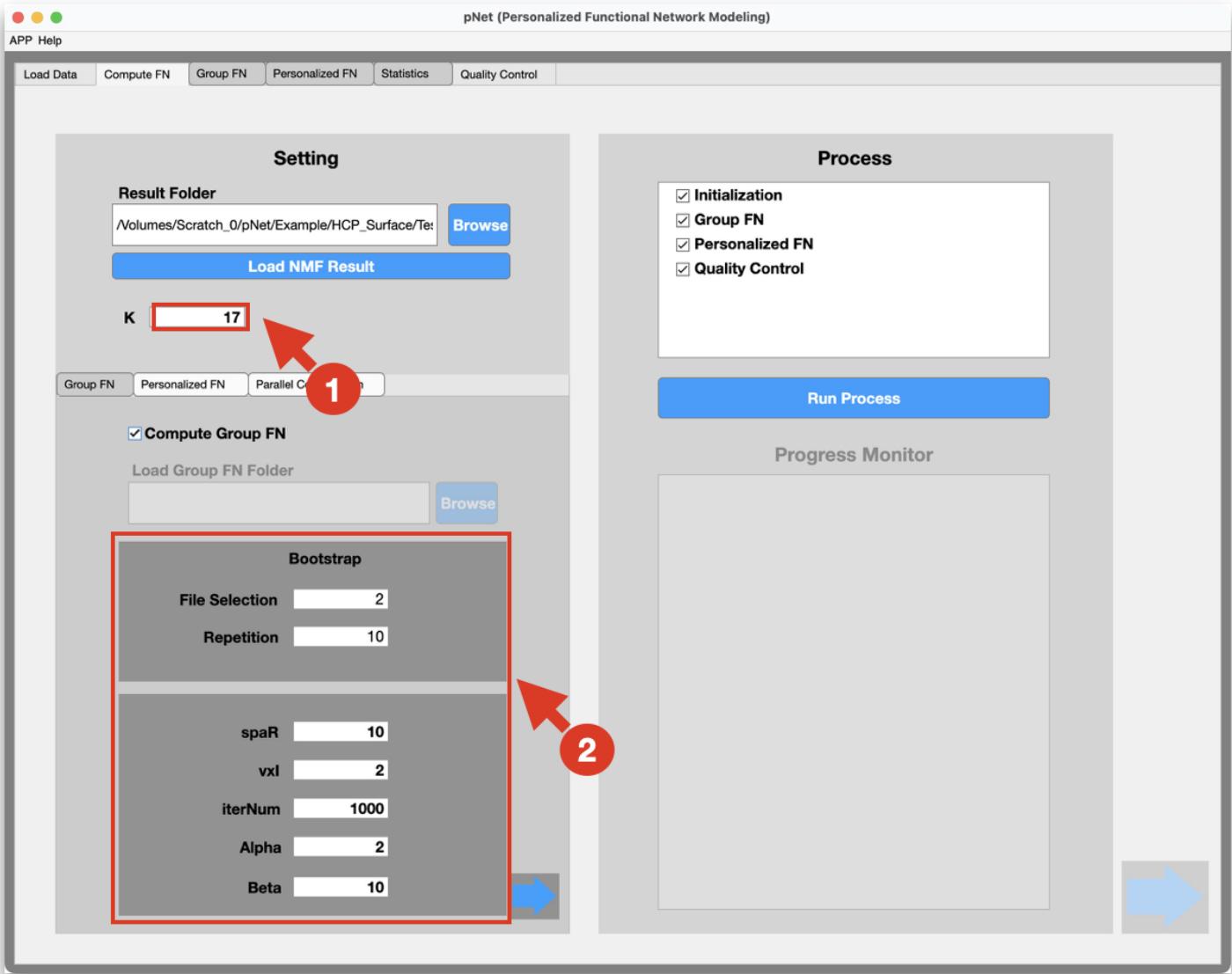

Setup the group-level FN computation parameters.

- (1) Set a desired number of FNs. 7 and 17 have been widely used to get large scale FNs.
- (2) Set additional parameters if preferred. The bootstrap section has two parameters that affect how the fMRI scans will be randomly sampled. The parameter "File Selection" is the number of fMRI scans will be selected as a subset to compute group-level FNs. The parameter "Repetition" is the number of runs to randomly sample a small subset and get multiple different versions group-level FN to fuse into a robust estimation of group-level FNs. Other parameters are for advanced users. Their details can be found in the manuscript.


2. Load a precomputed group-level FN file

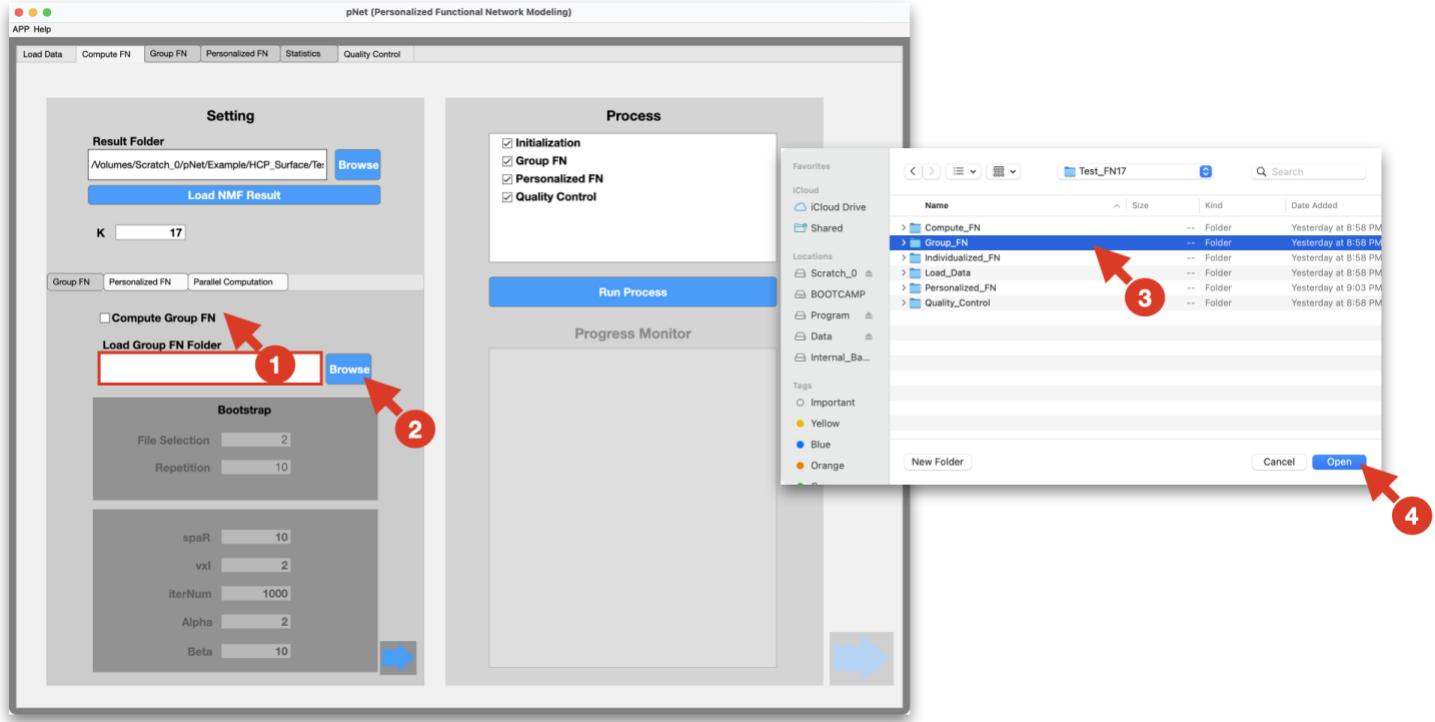

Load precomputed group-level FN results in a folder.

- (1) Click the check box “Compute Group FN” to disable the group-level FN computation.
- (2) Click the browse button or type in the directory of the folder containing group FN results in the input field (highlighted in the red box).
- (3-4) Select the group-level FN folder, and click open button to load the result.


Setup personalized FN computation

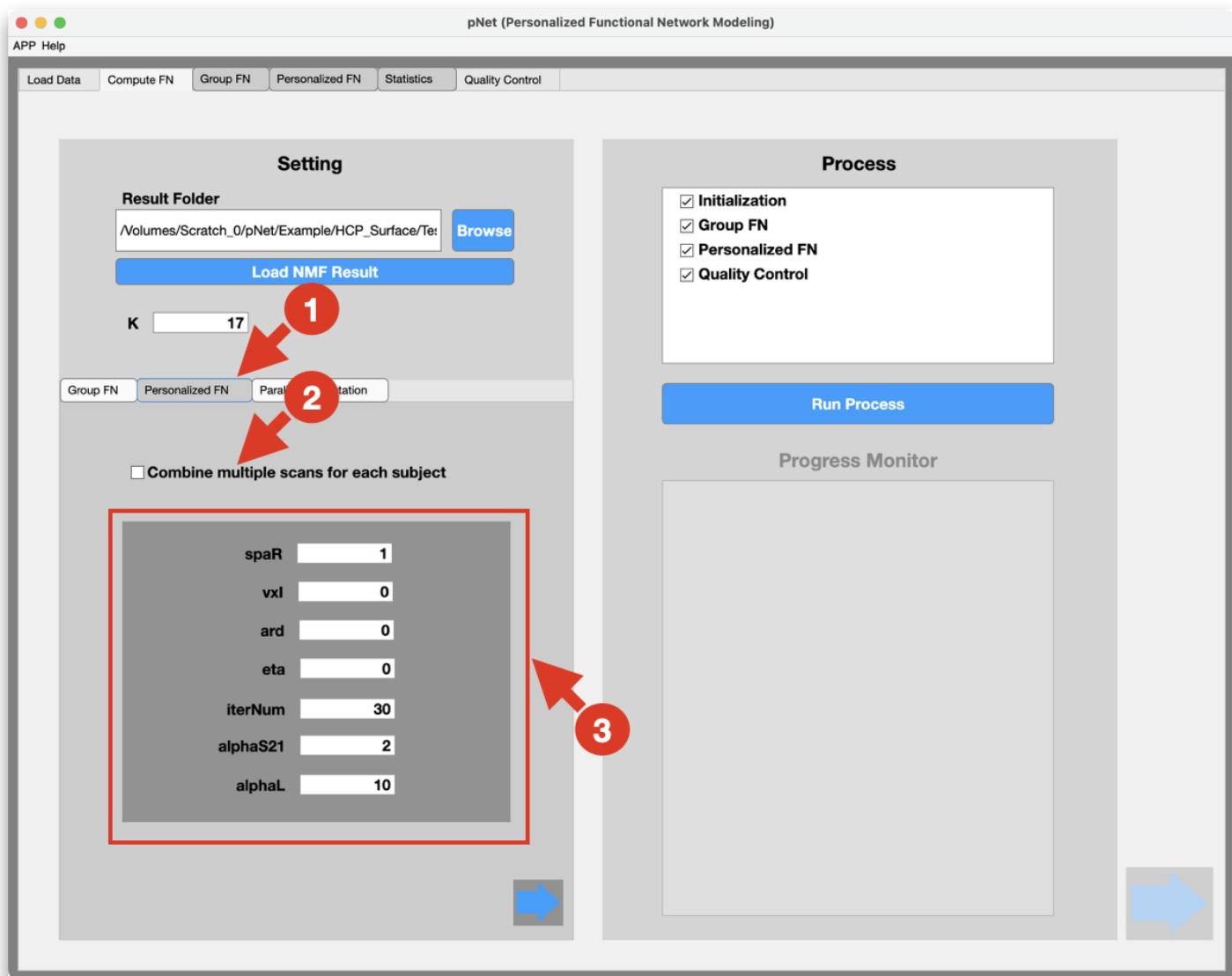

Setup parameters for personalized FN computation.

- (1) Navigate to the module “Personalized FN”.
- (2) Click the check box to combine multiple fMRI scans for each subject, if preferred.
- (3) Additional parameters for personalized FN computation, only intended for advanced users.

Setup parallel computation

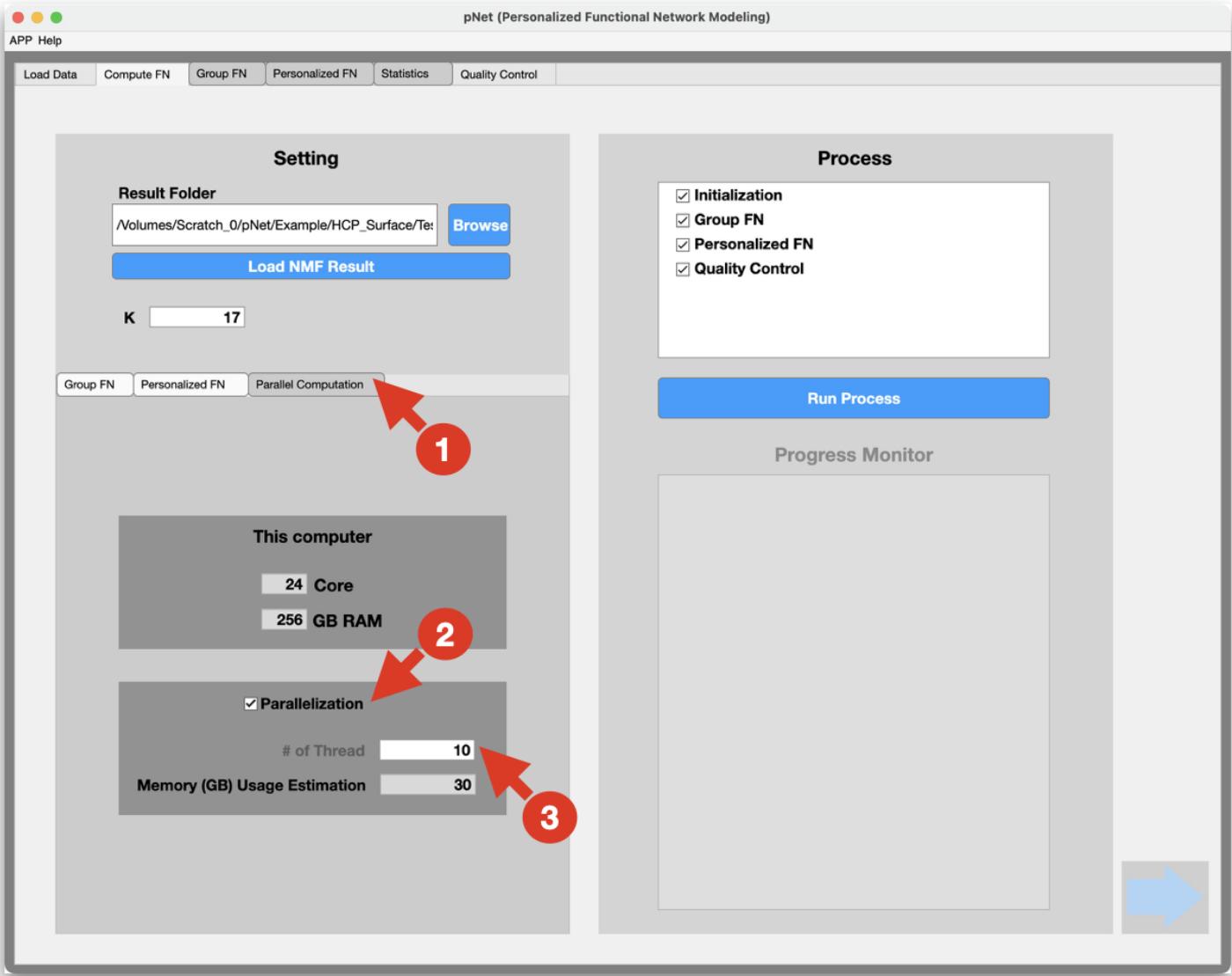


Setup parameters for parallel computation.

- (1) Navigate to the parallel computation module.
- (2) Enable parallelization by clicking the check box.
- (3) Set a desired number of threads for parallel computation. Please make sure that the estimated memory usage is within the hardware limit, which is displayed in the above section "This computer".

Run the FN computation

Once all the required settings are done, the button "Run Process" on the right of the "Compute FN" module will be enabled. Simply click it run all the computation and visualization for personalized FNs.

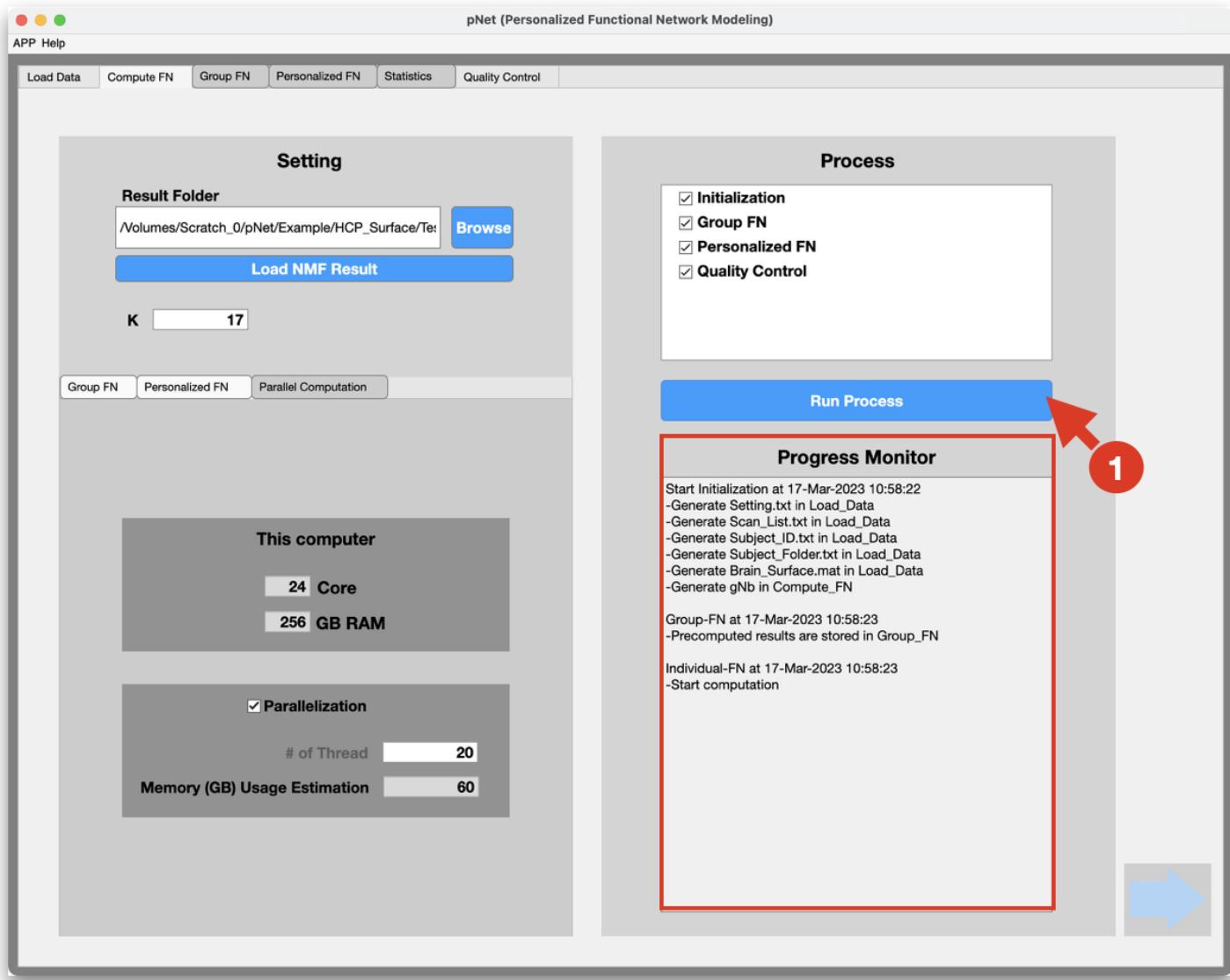

Run the personalized FN computation.


- (1) Click the button “Run Process”.
- (2) Progress monitor (highlighted in the red box) will show the progress of each major steps of initialization, group-level FN computation, personalized FN computation and visualization. The UI will automatically jump to the next module “Group FN” to show group-level FN results.

Check group-level FN

This toolbox provides precomputed figures of group-level FNs for fast navigation, and interactive display to check results according to users' preference.


Surface data

1. Precomputed visualization

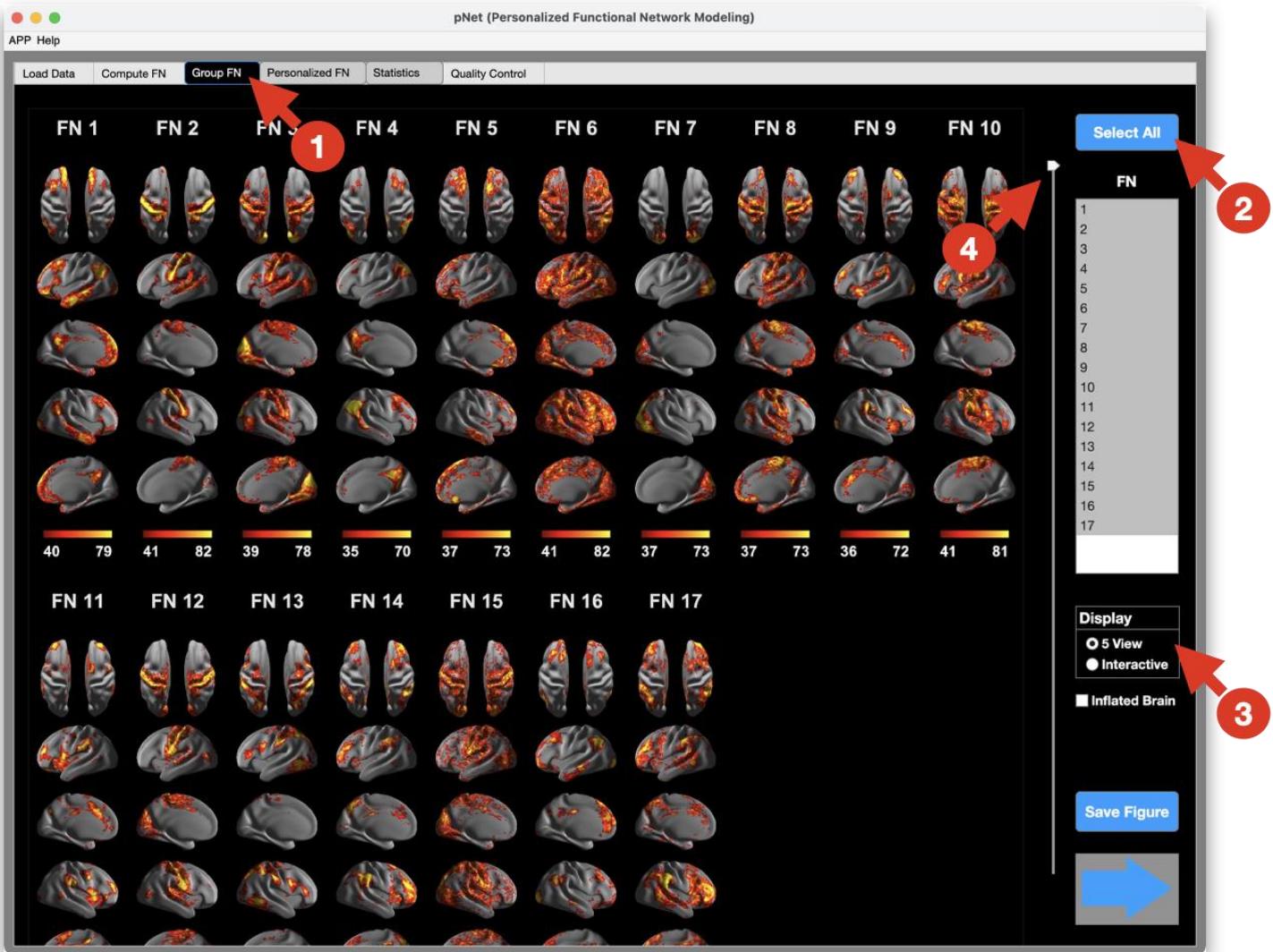

Check all group-level FNs.

- (1) Navigate to group FN module. In default, all group-level FNs will be displayed in one large figure.
- (2) Click the select all button, if needed.
- (3) Select the 5 view display option, if needed.
- (4) Drag the scroll bar to navigate through the entire figure.

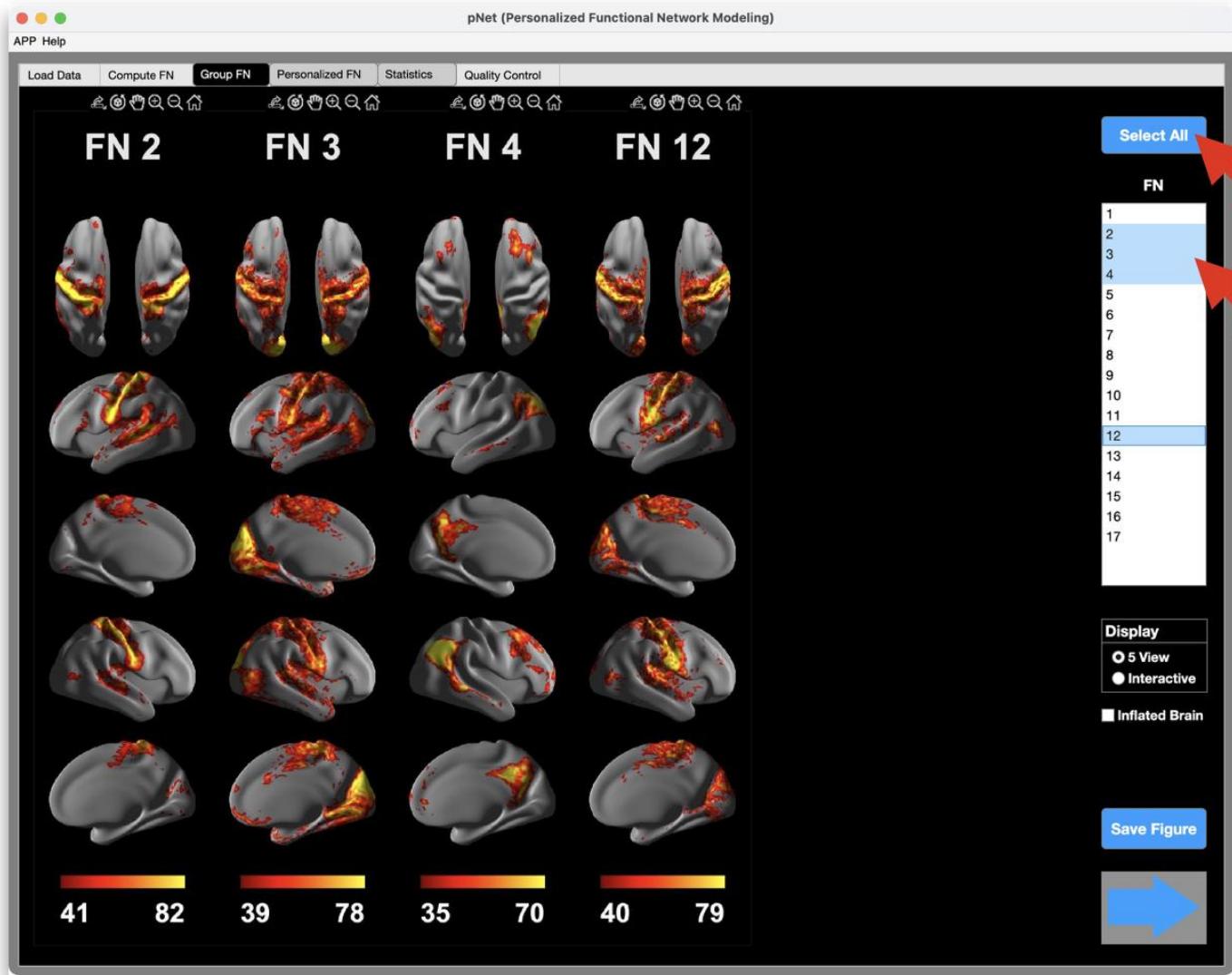


Check selected group-level FNs.

- (1) Select one or multiple group-level FNs. Hold command key in macOS, or hold control key in Windows or Linux to select multiple group-level FNs.
- (2) Click select all button to go back to the default view of all group-level FNs selected.


2. Interactive visualization

This interactive visualization allows the user check the group-level FN from different view angles.

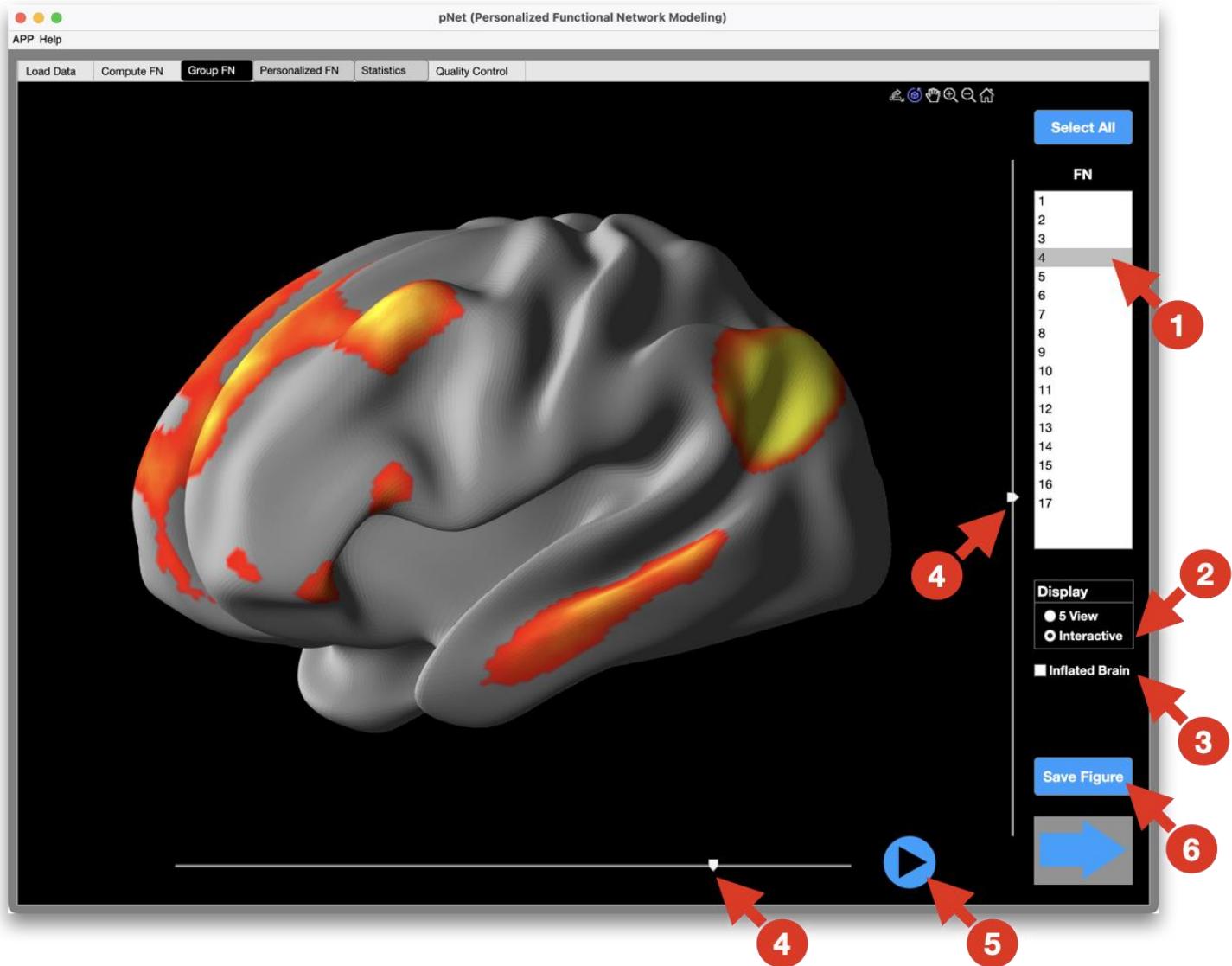

Check a single group-level FN on the brain surface.

- (1) Select one group-level FN.
- (2) Select the interactive display option.
- (3) Click the check box to overlay results on an inflated brain surface, if the inflated brain surface is available.
- (4) Drag the vertical or horizontal scroll bar to change the view angle.
- (5) Start automatic change of the horizontal scroll bar to have a constantly rotating view of the brain.
- (6) Save this figure.

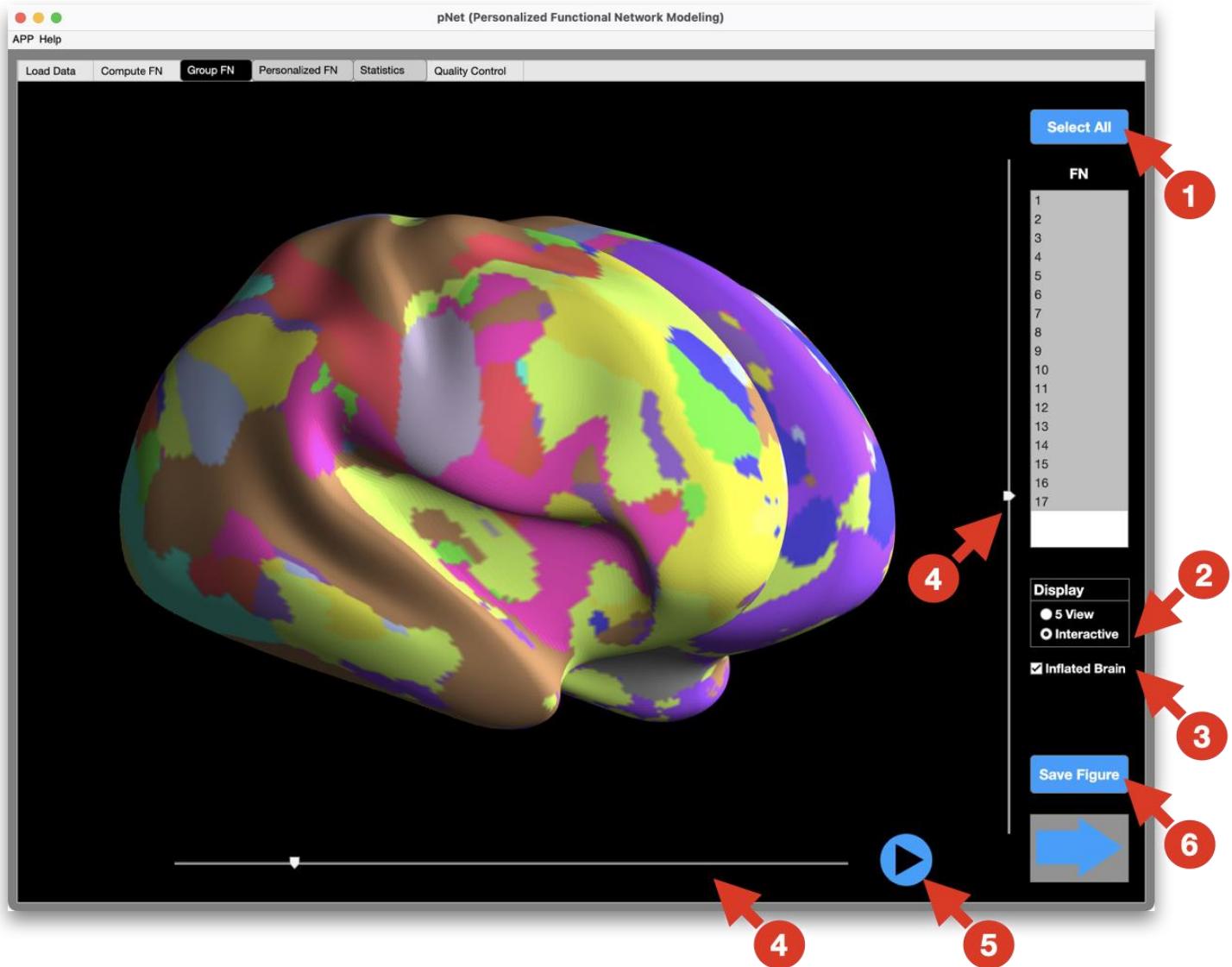


Check the functional atlas for binarized group-level FNs.

- (1) Click the select all button to include all group-level FNs.
- (2) Select the interactive display option.
- (3) Click the check box to overlay results on an inflated brain surface, if the inflated brain surface is available.
- (4) Drag the vertical or horizontal scroll bar to change the view angle.
- (5) Start automatic change of the horizontal scroll bar to have a constantly rotating view of the brain.
- (6) Save this figure.


Volume data

1. Precomputed visualization

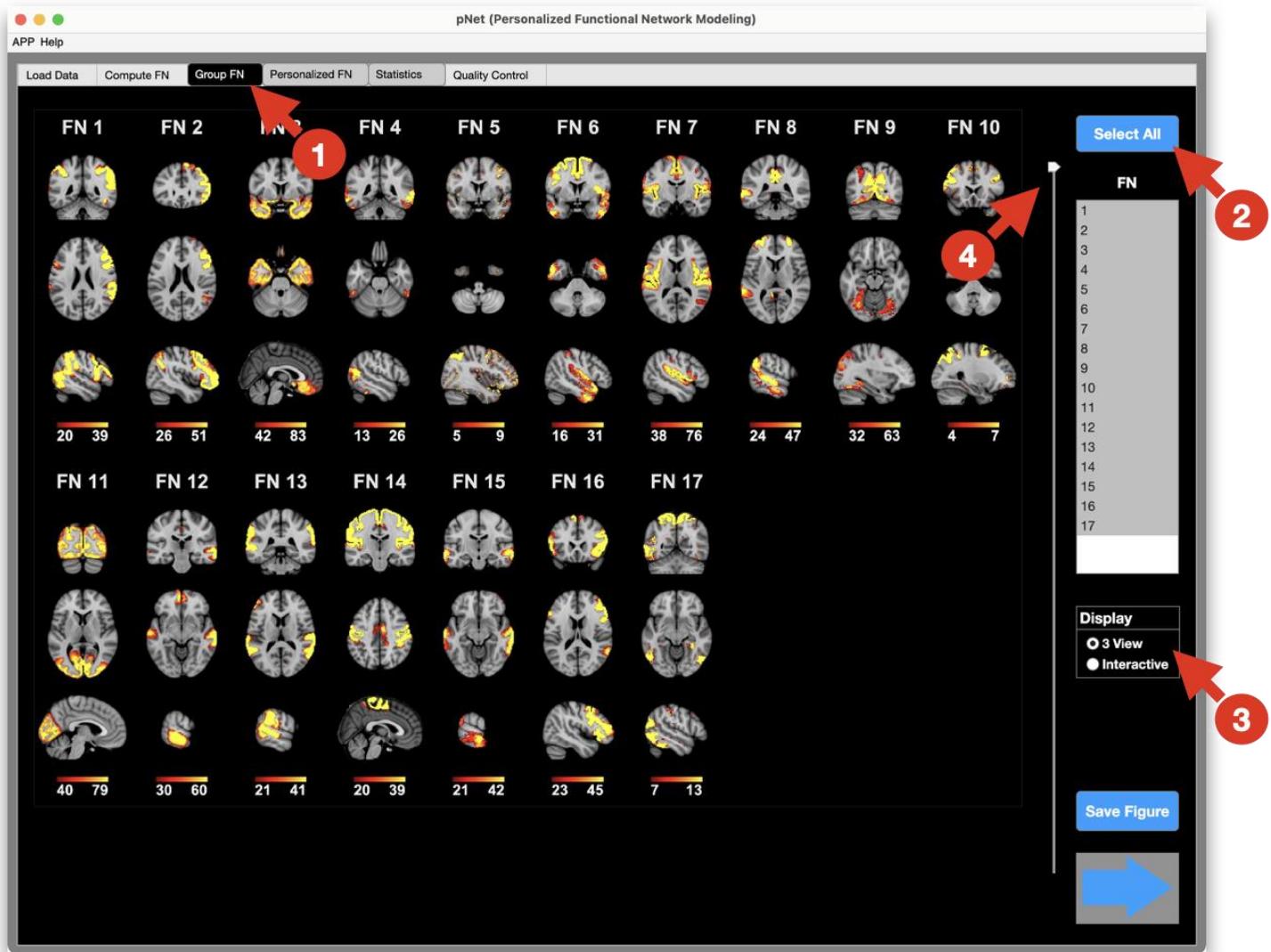

Check all group-level FNs.


- (1) Navigate to group FN module. In default, all group-level FNs will be displayed in one large figure.
- (2) Click the select all button, if needed.
- (3) Select the 3 view display option, if needed.
- (4) Drag the scroll bar to navigate through the entire figure.

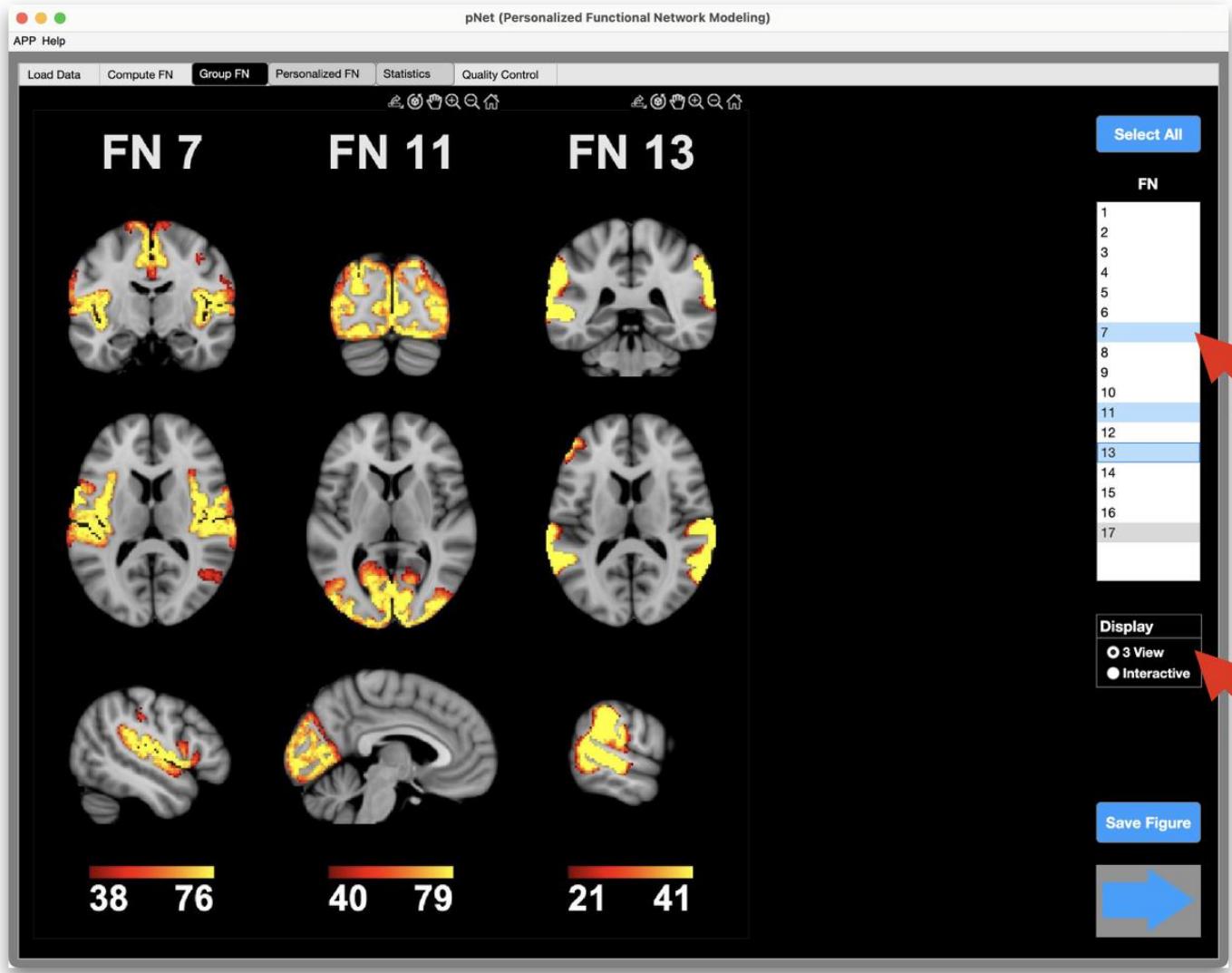

Check selected group-level FNs.


- (1) Select one or multiple group-level FNs. Hold command key in macOS, or hold control key in Windows or Linux to select multiple group-level FNs.
- (2) Click select all button to go back to the default view of all group-level FNs selected.

2. Interactive visualization

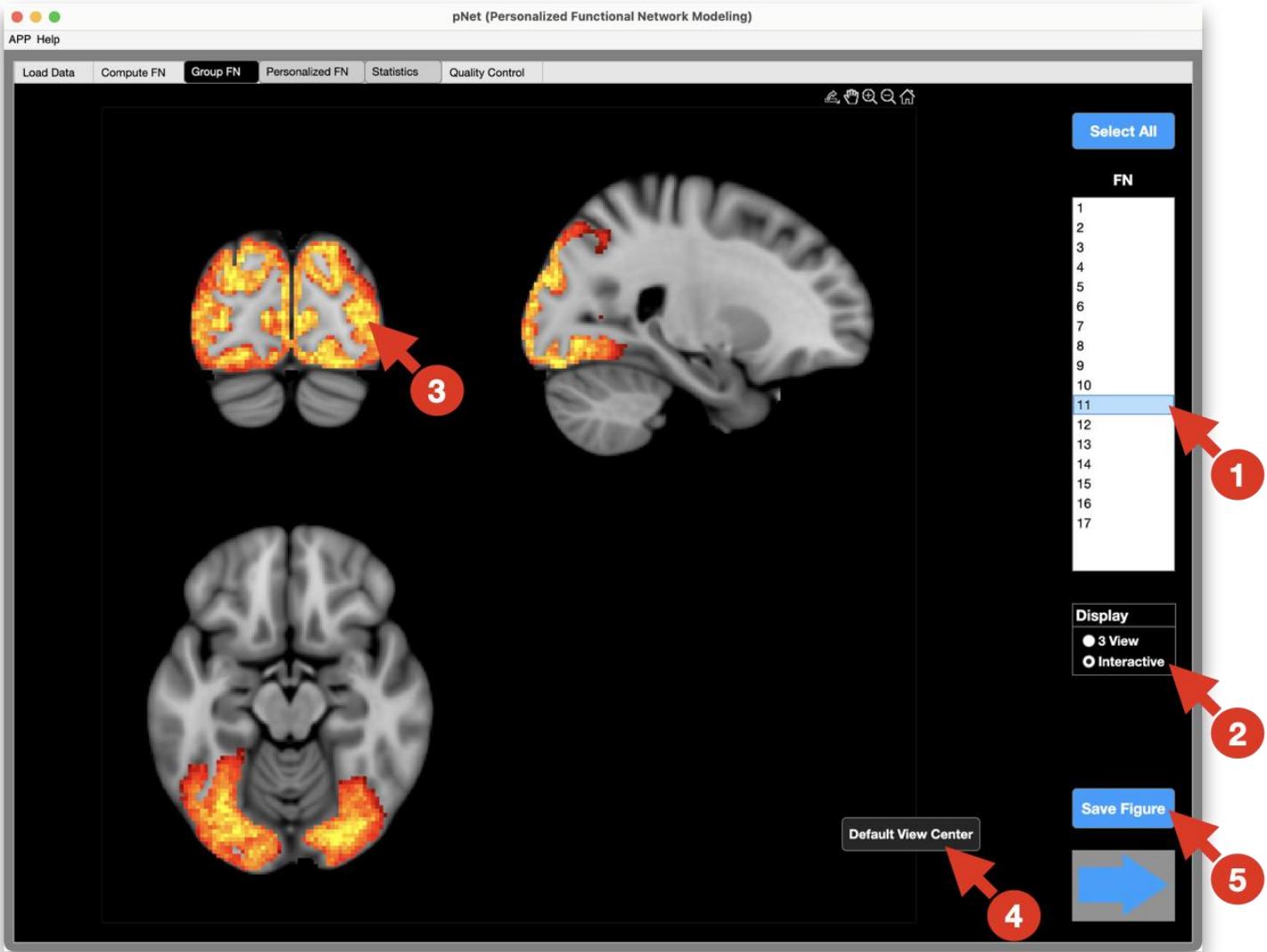

Check a single group-level FN

- (1) Select one group-level FN.
- (2) Select the interactive display option.
- (3) Click on the display to change the view center.
- (4) Click the default view center button to go back to the default view center.
- (5) Save this figure.

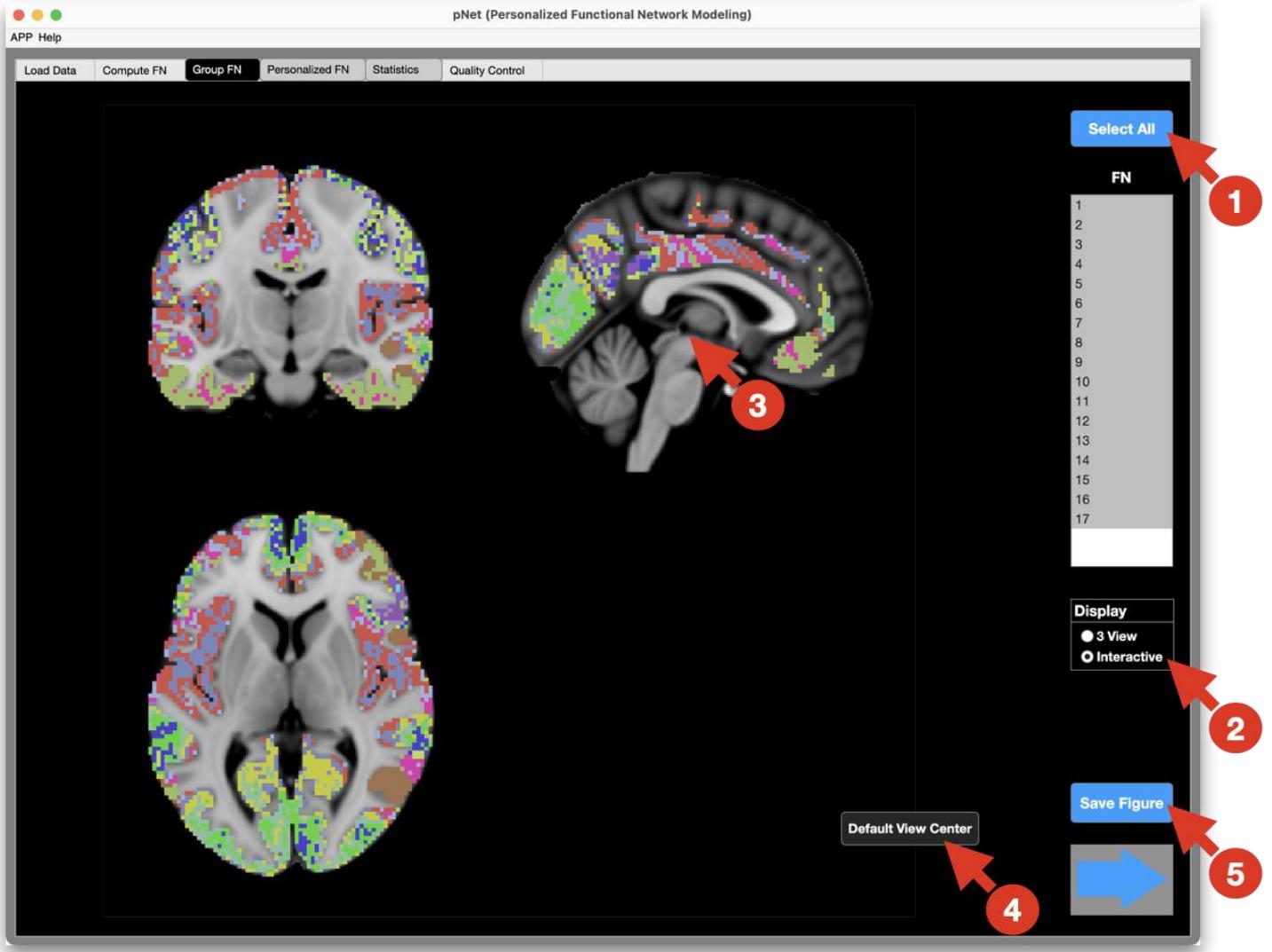

Check functional atlas using binarized group-level FN results


- (1) Click the select all button to include all group-level FNs.
- (2) Select the interactive display option.
- (3) Click on the display to change the view center.
- (4) Click the default view center button to go back to the default view center.
- (5) Save this figure.

Check personalized FN

This toolbox provides precomputed figures of personalized FNs for fast navigation, and interactive display to check results according to users' preference.


Surface data

1. Precomputed visualization

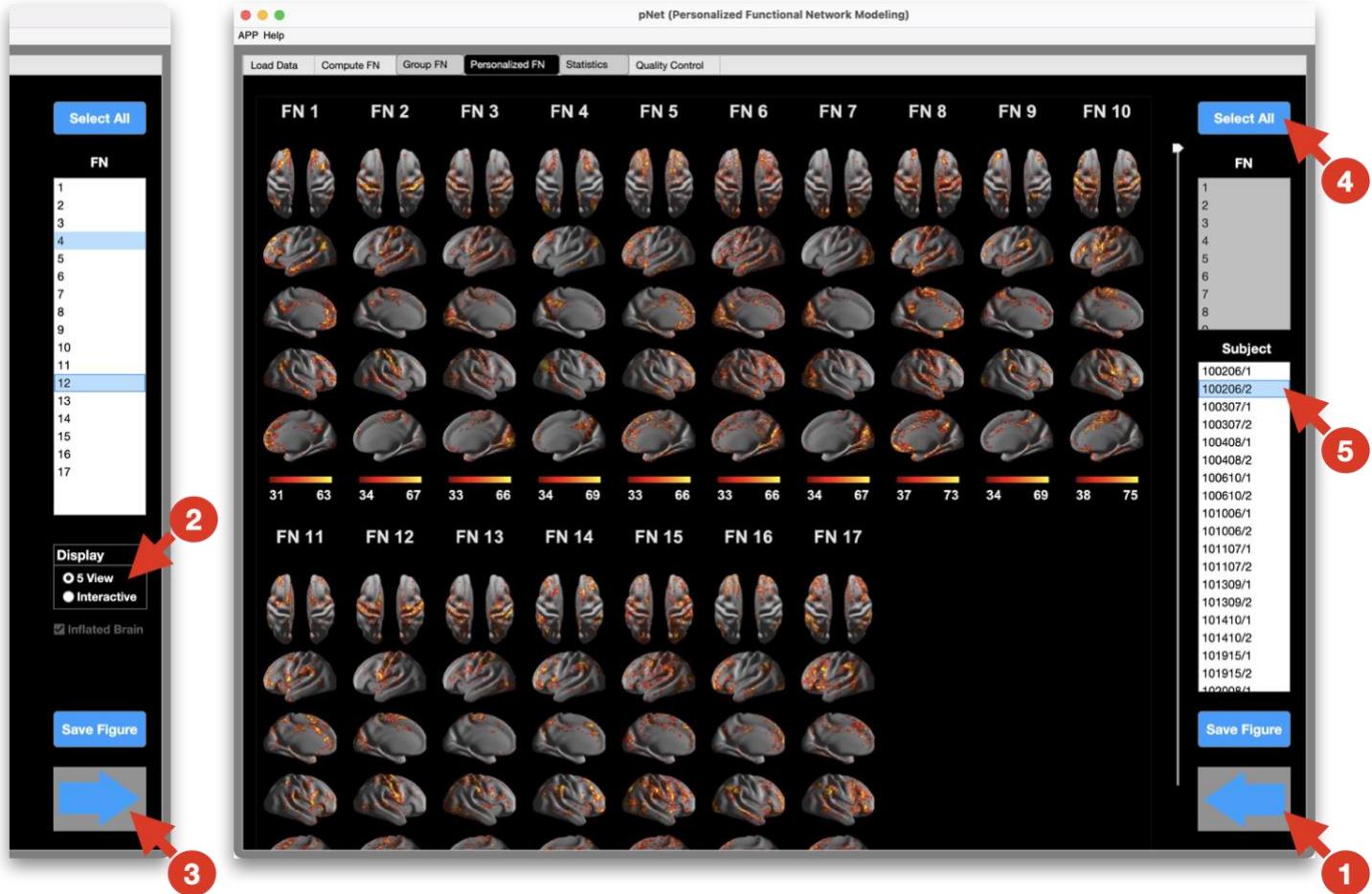

Check precomputed figures for personalized FNs.


- (1) Navigate to the group FN module.
- (2) Choose the 5 view display mode.
- (3) Navigate to the personalized FN module.
- (4) Click the select all button.
- (5) Select another subject. Users can also use the up or down arrow key to select the subject above or below.

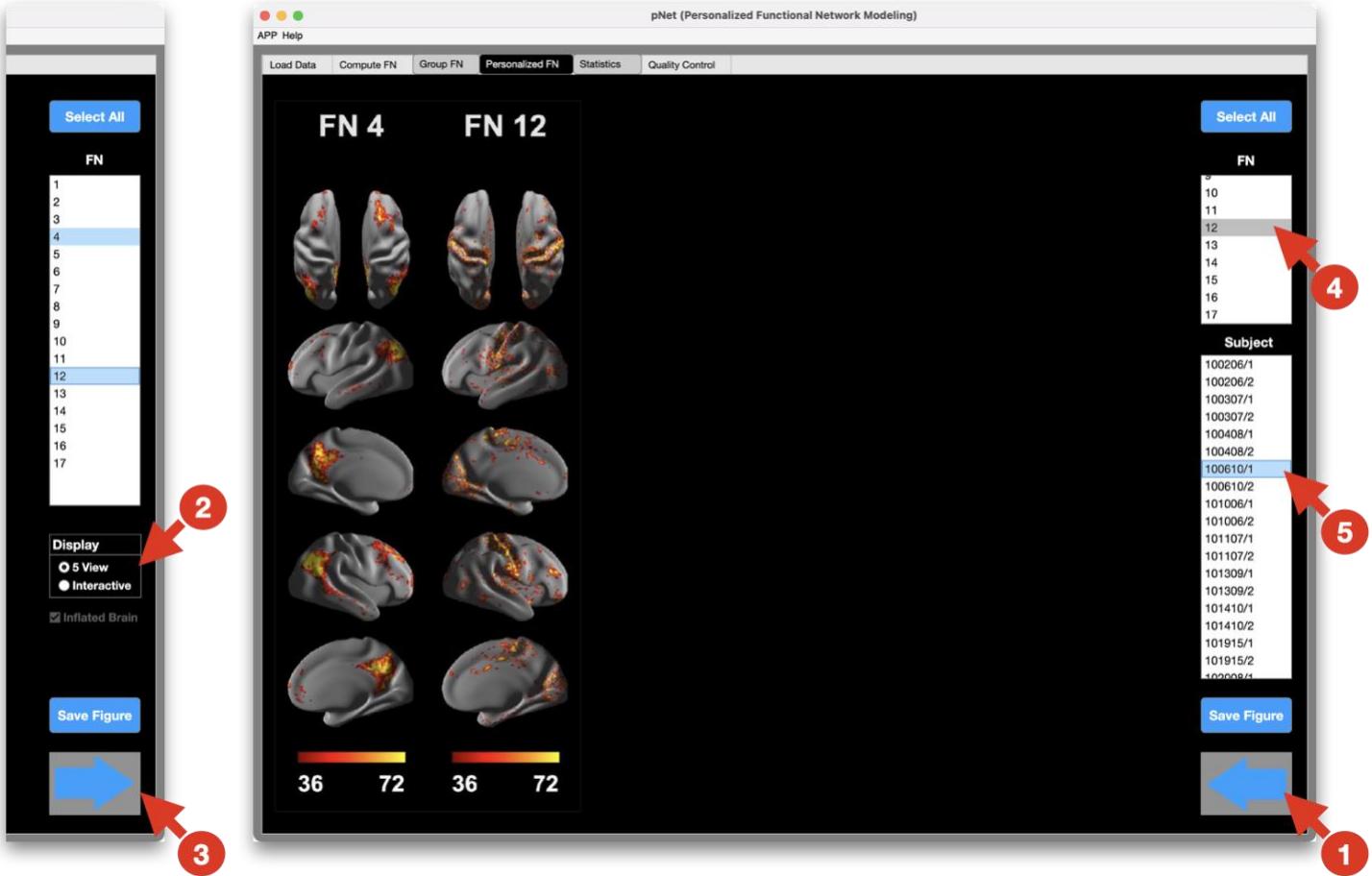

Check selected personalized FNs.


- (1) Navigate to the group FN module
- (2) Select the 5 view display option.
- (3) Navigate to the personalized FN module.
- (4) Select one or multiple personalized FNs. Hold command key in macOS, or hold control key in Windows or Linux to select multiple personalized FNs.
- (5) Select another subject. Users can also use the up or down arrow key to select the subject above or below.

2. Interactive visualization

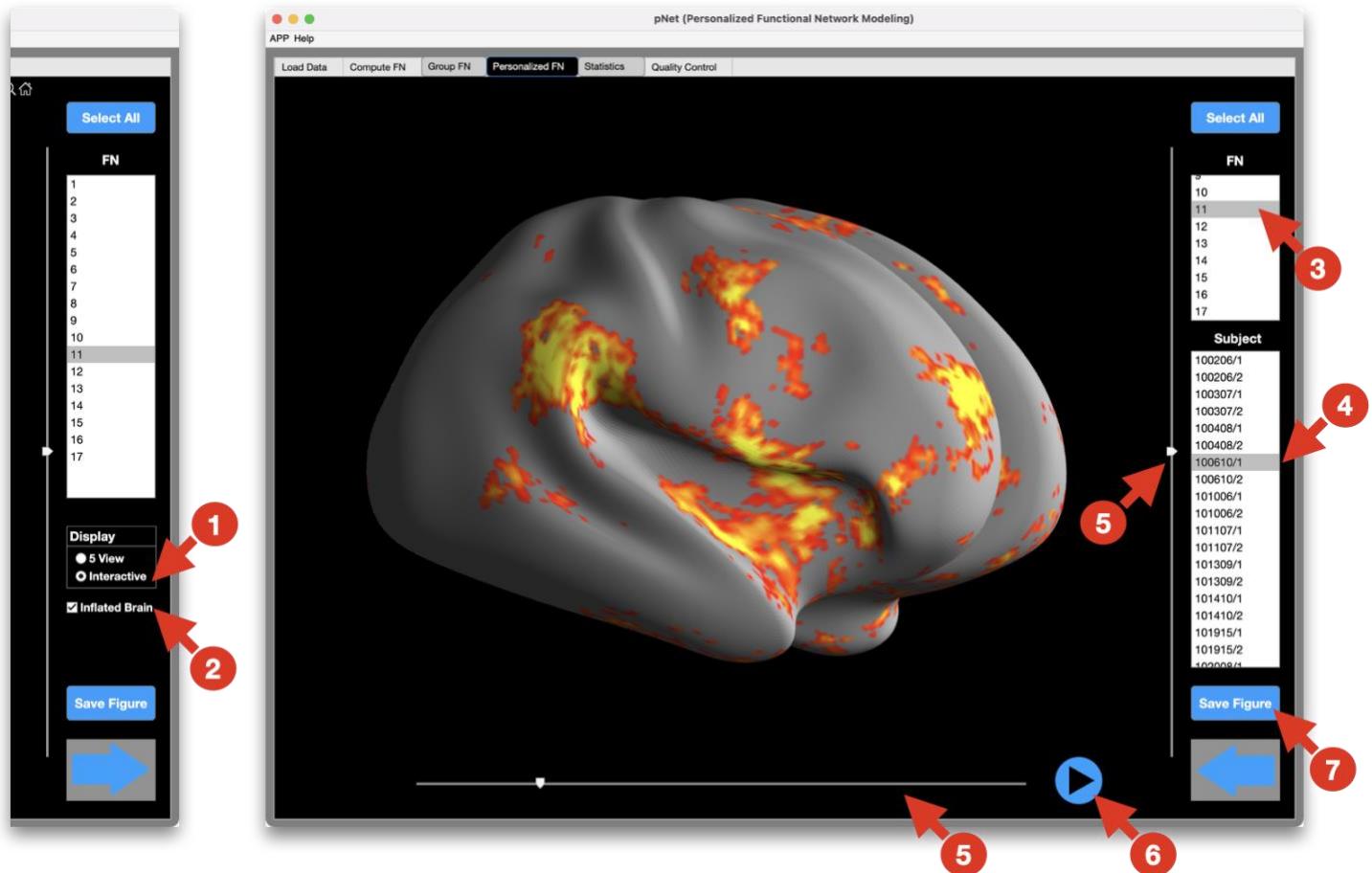

Check a single personalized FN on the brain surface.

- (1) Navigate to the group FN module, and select the interactive display option.
- (2) Select the inflated brain if desired.
- (3) Navigate to the personalized FN module, select one desired personalized FN for display.
- (4) Select another subject. Users can also use the up or down arrow key to select the subject above or below.
- (5) Drag the vertical or horizontal scroll bar to change the view angle.
- (6) Start automatic change of the horizontal scroll bar to have a constantly rotating view of the brain.
- (7) Save this figure.

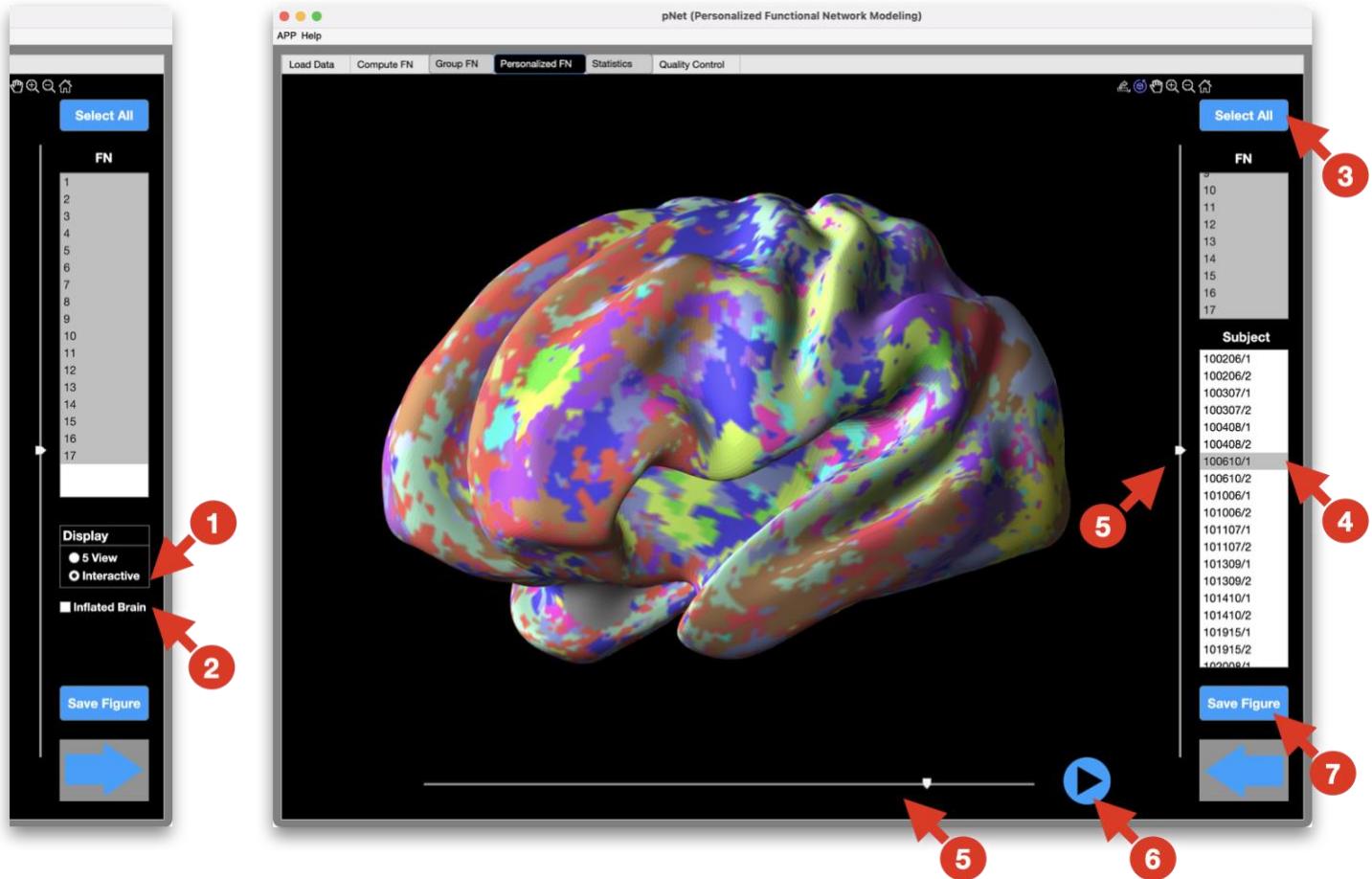


Check the functional atlas of a single subject.

- (1) Navigate to the group FN module, and select the interactive display option.
- (2) Select the inflated brain if desired.
- (3) Navigate to the personalized FN module, click the select all button if needed.
- (4) Select another subject. Users can also use the up or down arrow key to select the subject above or below.
- (5) Drag the vertical or horizontal scroll bar to change the view angle.
- (6) Start automatic change of the horizontal scroll bar to have a constantly rotating view of the brain.
- (7) Save this figure.


Volume data

1. Precomputed visualization

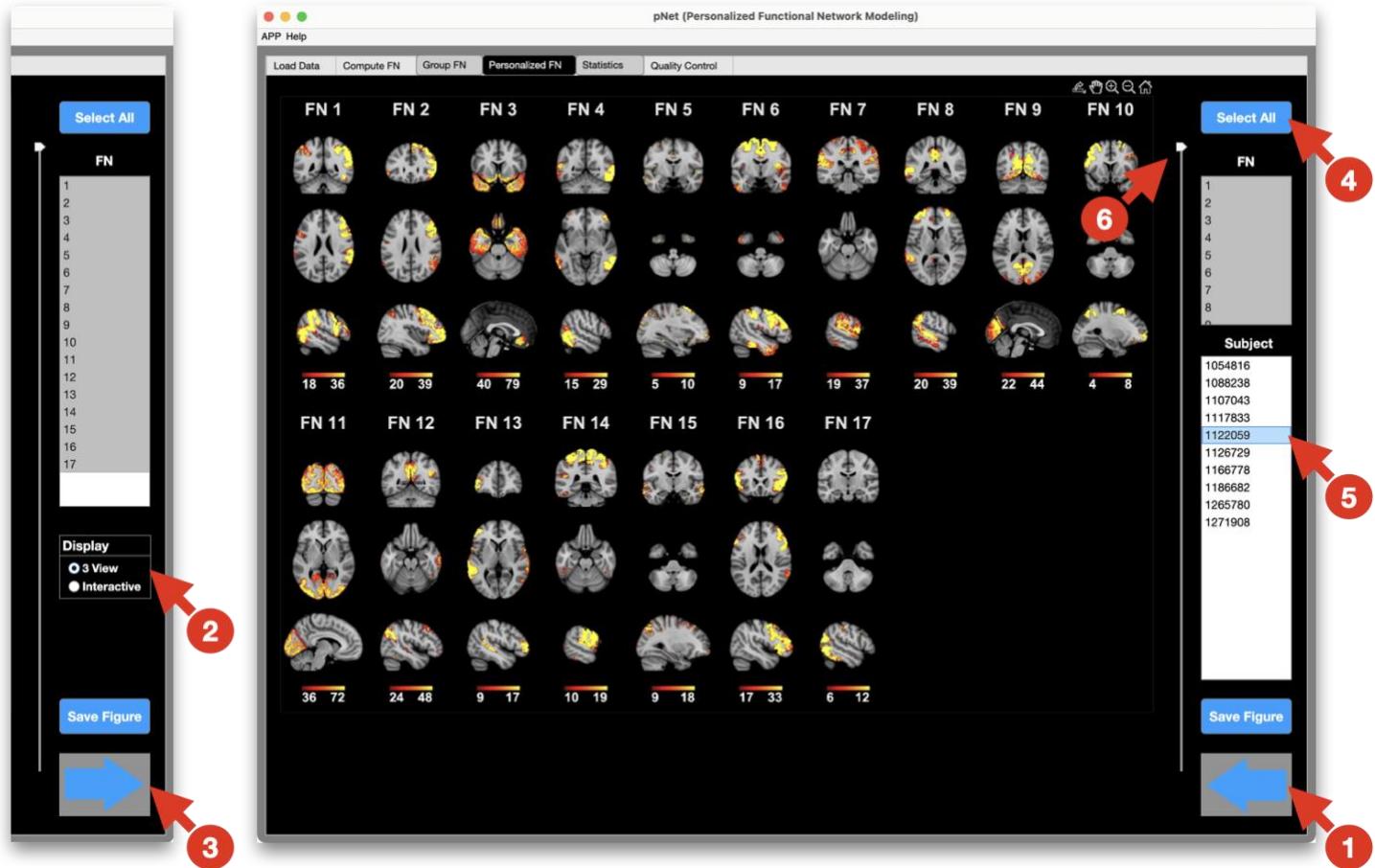

Check precomputed figures for personalized FNs.


- (1) Navigate to the group FN module.
- (2) Choose the 3 view display mode.
- (3) Navigate to the personalized FN module.
- (3) Click the select all button, if needed.
- (5) Select another subject. Users can also use the up or down arrow key to select the subject above or below.
- (6) Drag the vertical scroll bar, if needed.

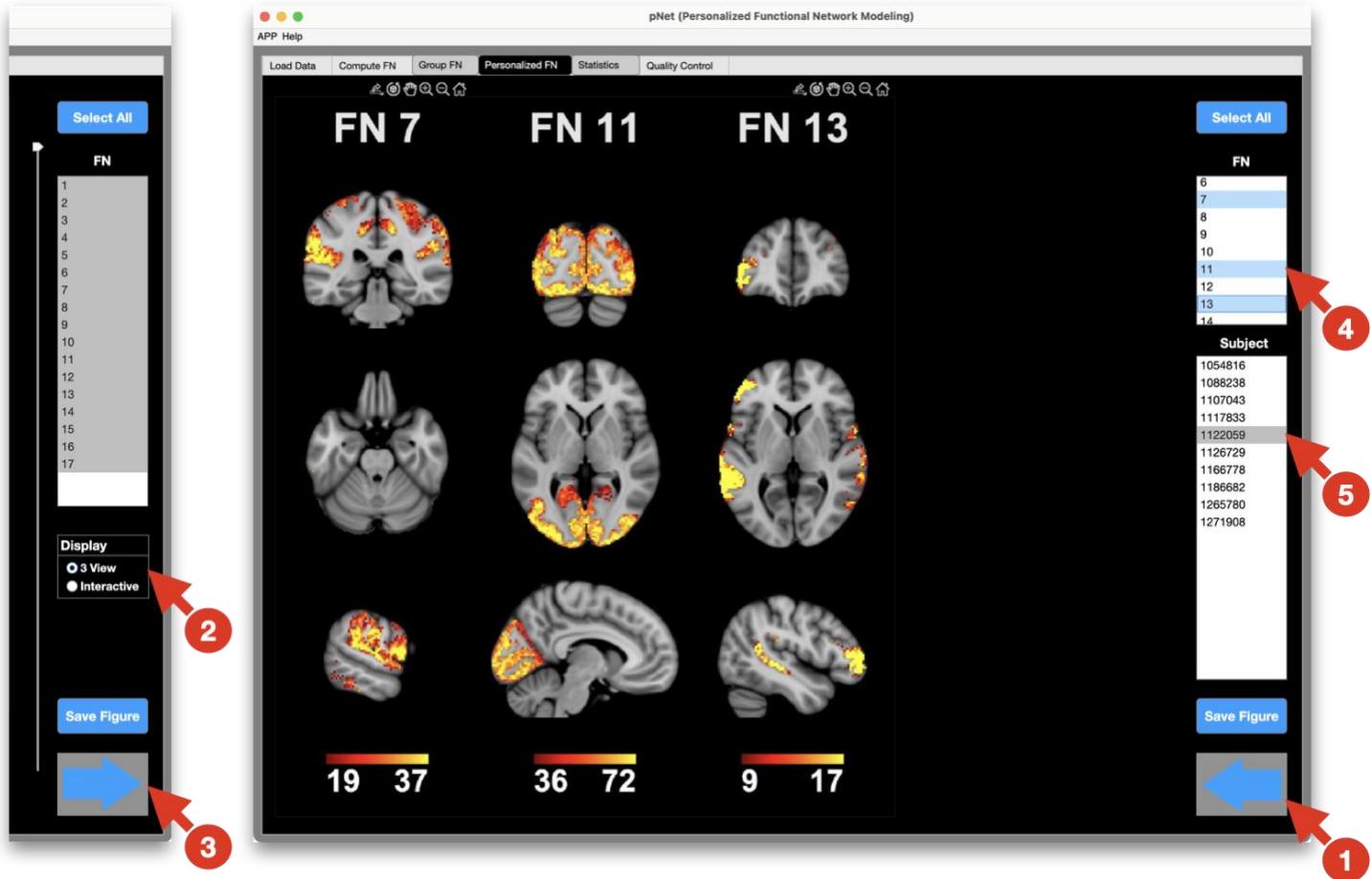

Check selected personalized FNs.


- (1) Navigate to the group FN module.
- (2) Select the 3 view display option.
- (3) Navigate to the personalized FN module.
- (4) Select one or multiple personalized FNs. Hold command key in macOS, or hold control key in Windows or Linux to select multiple personalized FNs.
- (5) Select another subject. Users can also use the up or down arrow key to select the subject above or below.

2. Interactive visualization

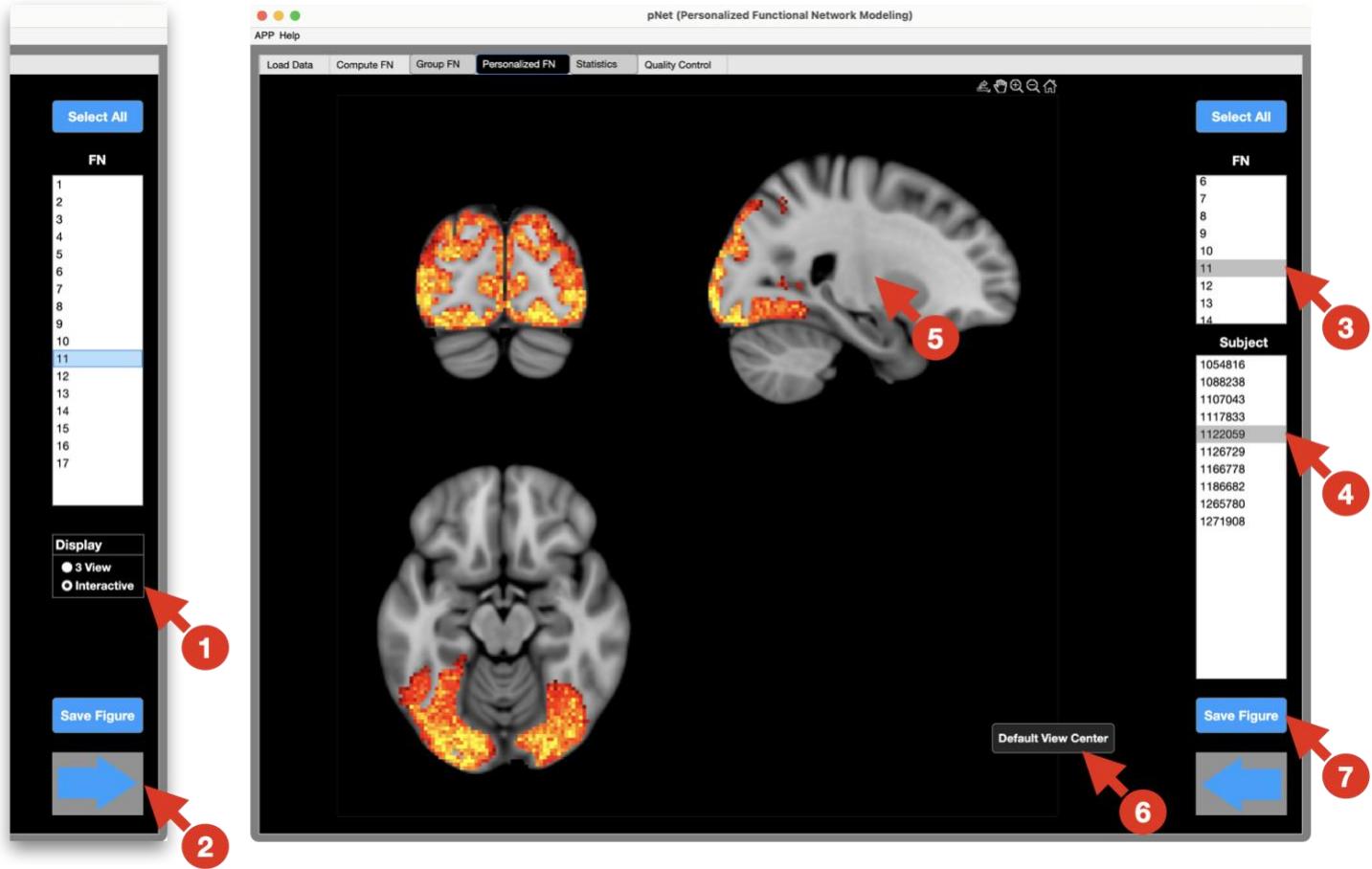

Check a single personalized FN.

- (1) Navigate to the group FN module, and select the interactive display option.
- (2) Navigate to the personalized FN module
- (3) Select one desired personalized FN for display.
- (4) Select another subject. Users can also use the up or down arrow key to select the subject above or below.
- (5) Click on the display to change the view center.
- (6) Click the default view center button to go back to the default view center.
- (7) Save this figure.

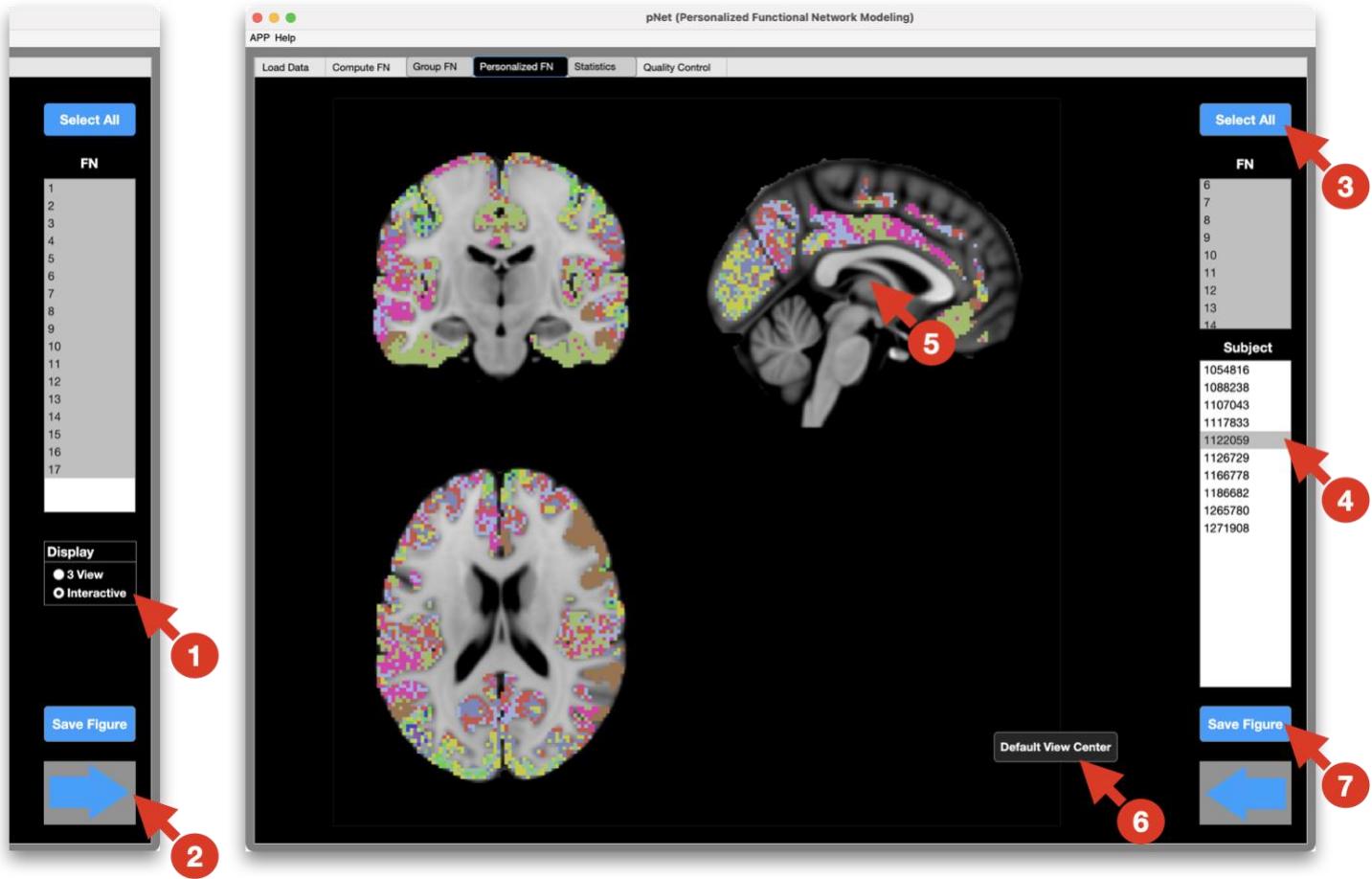

Check the functional atlas of a single subject.

- (1) Navigate to the group FN module, and select the interactive display option.
- (2) Navigate to the personalized FN module.
- (3) Click select all button, if needed.
- (4) Select another subject. Users can also use the up or down arrow key to select the subject above or below.
- (5) Click on the display to change the view center.
- (6) Click the default view center button to go back to the default view center.
- (7) Save this figure.

Compare group-level and personalized FN

Users can easily compare group-level and personalized FNs using the fast navigation buttons on the right bottom of the group-level and personalized FN modules. The display settings for group-level and personalized FNs are synchronized automatically for convenience.




Fast navigation between group-level and personalized FN modules.

- (1) The right arrow button on the group-level FN module navigates to the personalized FN module.
- (2) The left arrow button on the personalized FN module navigates to the group-level FN module.

Statistical Analysis

This toolbox has implemented one-sample t-test to show the significance of each personalized FN, and two-sample t-test to show the significant difference of each personalized FN in two groups. The visualization options are as same as those for group-level and personalized FNs. In addition, the visualization for statistical results are synchronized with that for group-level and personalized FNs. It allows users to navigate to the group-level and personalized FN modules for fast comparisons.

One-sample t-test

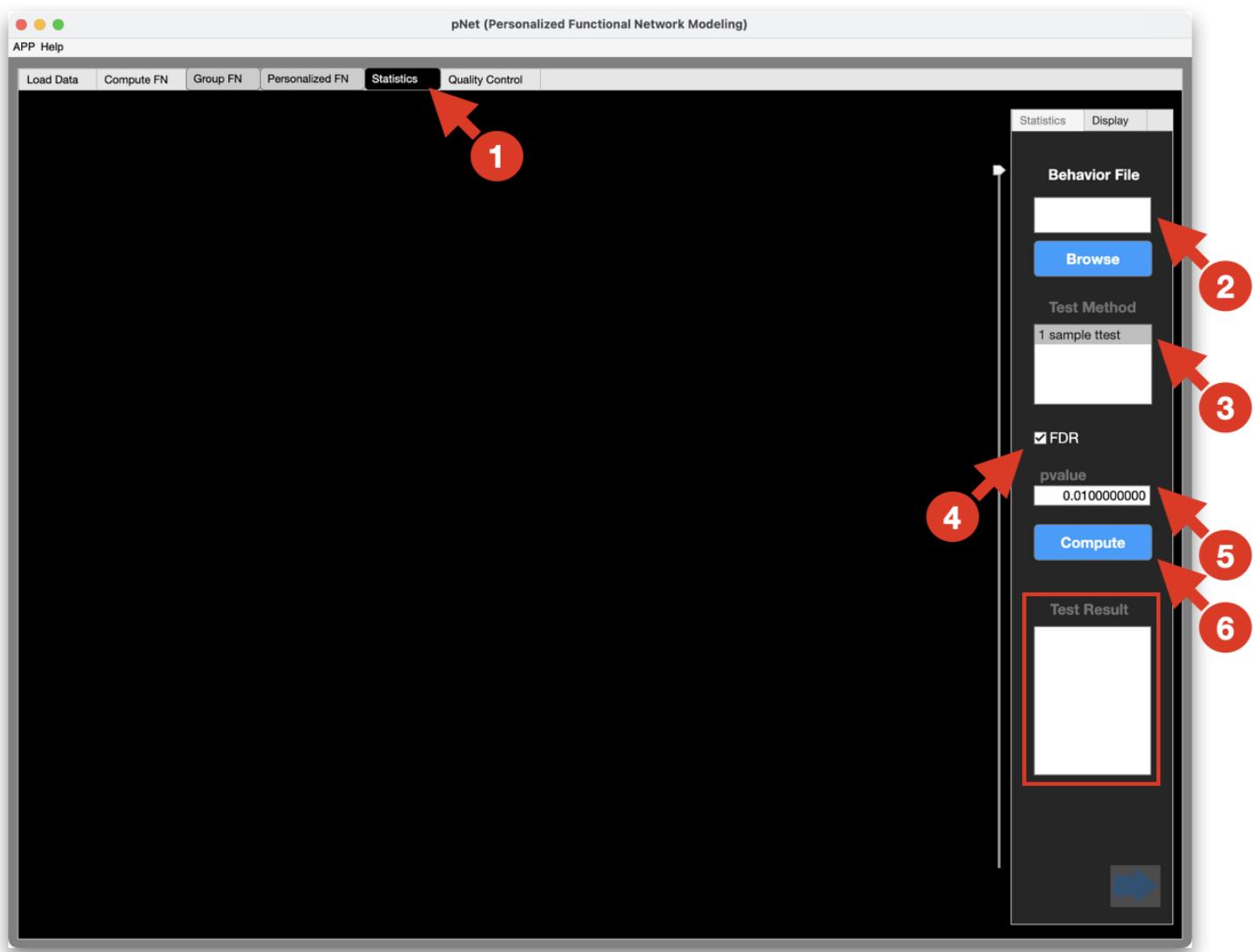


Setup one-sample t-test.

- (1) Navigate to the statistics module.
- (2) Make sure that the behavior file is empty.
- (3) Select the 1 sample t-test method.


- (4) Check FDR correction, if needed.
- (5) Set the p-value threshold.
- (6) Click the compute button to run the statistical analysis. A popup window will show the progress of statistical computation, and another window will show the progress of result visualization.

Once the computation and visualization finish, a new test result will be added into the test result list box.

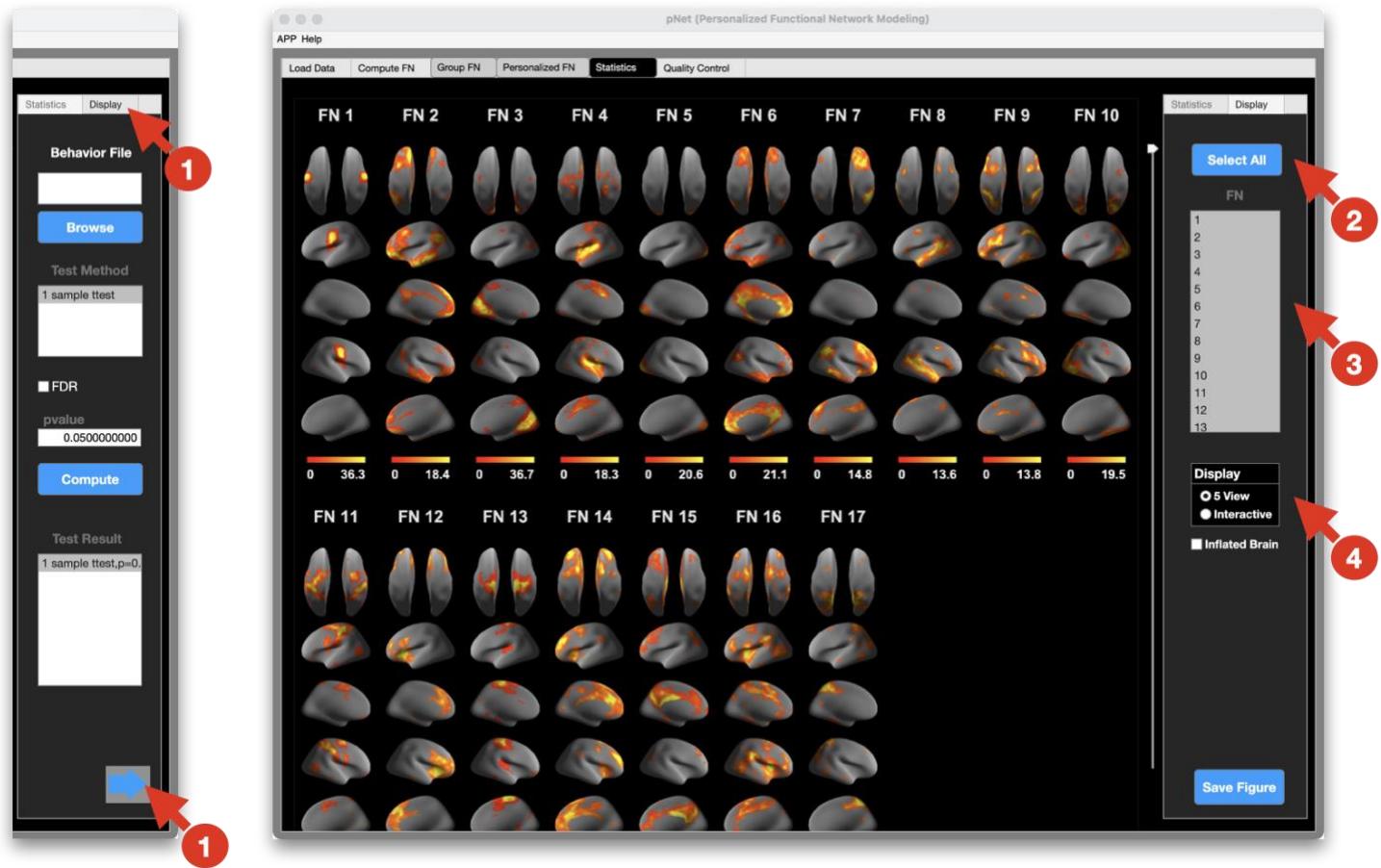

Check precomputed visualization.


- (1) Navigate to the display module.
- (2) Click the select all button, if needed.
- (3) Select one or multiple FNs.
- (4) Select the 5 view display option, if needed.

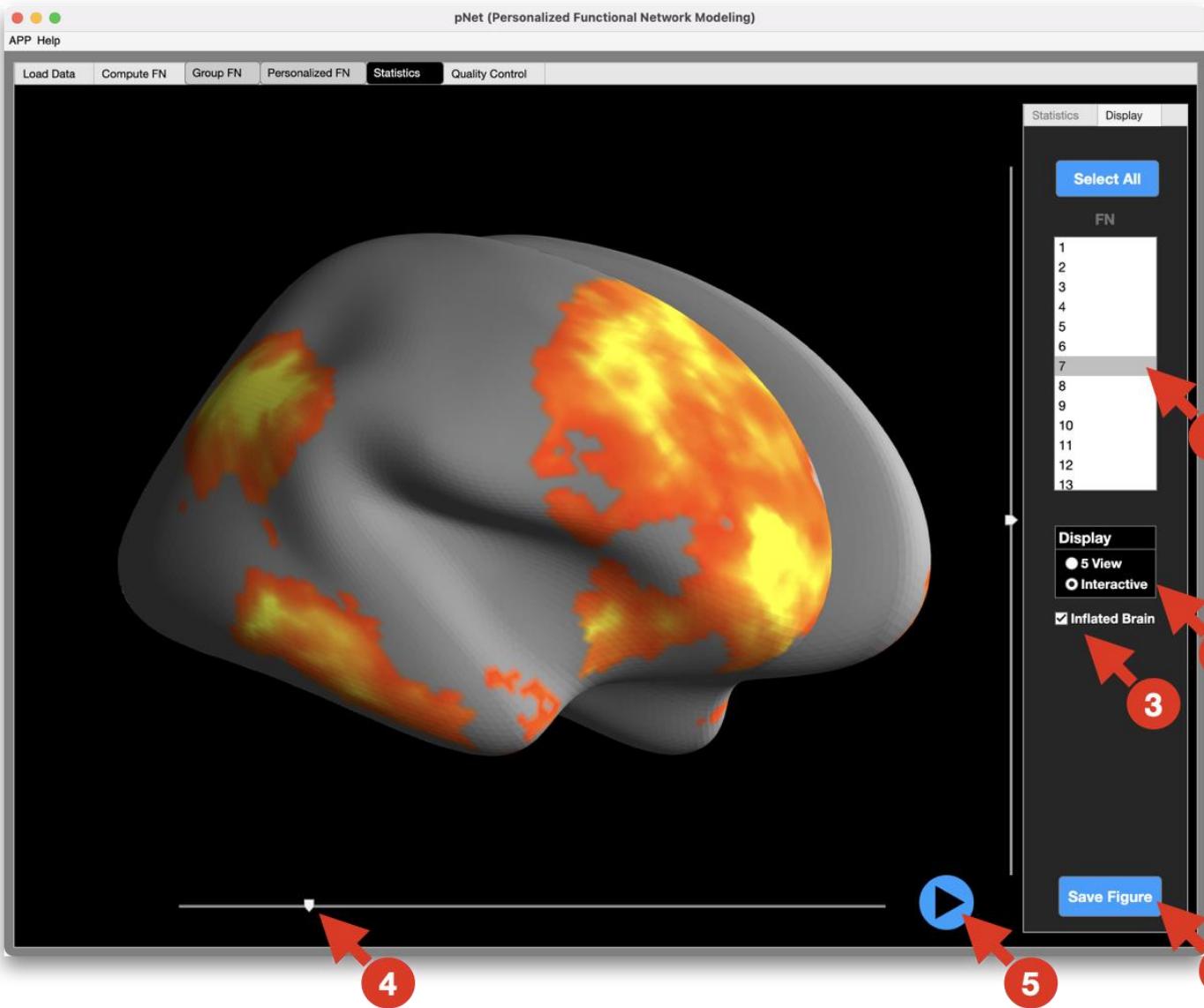

Check interactive visualization.

- (1) Select one FN.
- (2) Select the interactive display option.
- (3) Choose the inflated brain, if preferred.
- (4) Drag the scroll bar to change the view angle.
- (5) Click the play button to have a constantly rotating view of the brain.
- (6) Click the save figure button to save this figure.


Two-sample t-test

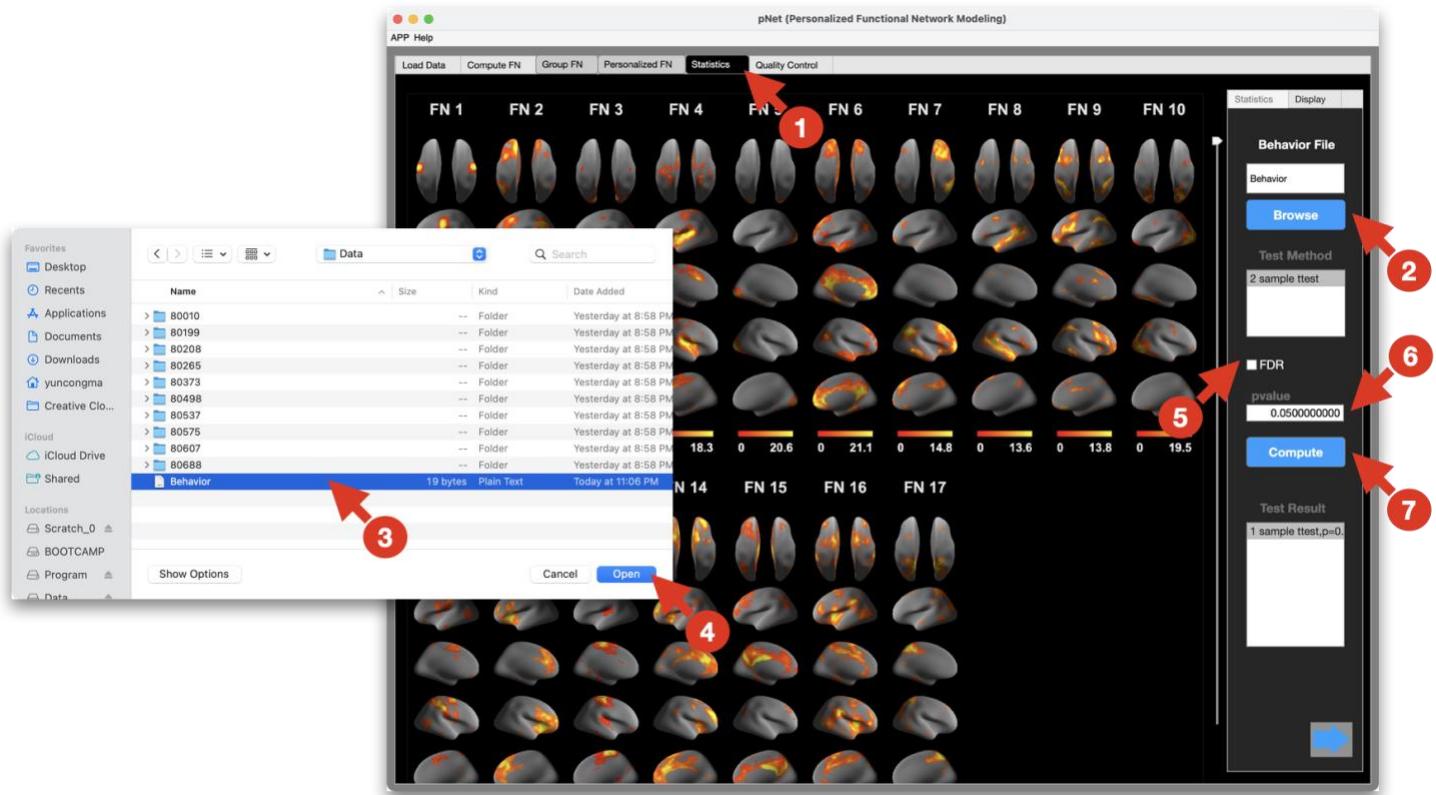

Setup the two-sample t-test

- (1) Navigate to the statistics module.
- (2) Click the browse button to read a behavior file.
- (3) Select the behavior file in txt format, which contains a column of 0-1 to specify the two groups.
- (4) Click open button.
- (5) Select the FDR correction, if preferred.
- (6) Setup the p-value threshold.
- (7) Click the compute button to run the statistical analysis. A popup window will show the progress of statistical computation, and another window will show the progress of result visualization.

Once the computation and visualization finish, a new test result will be added into the test result list box.

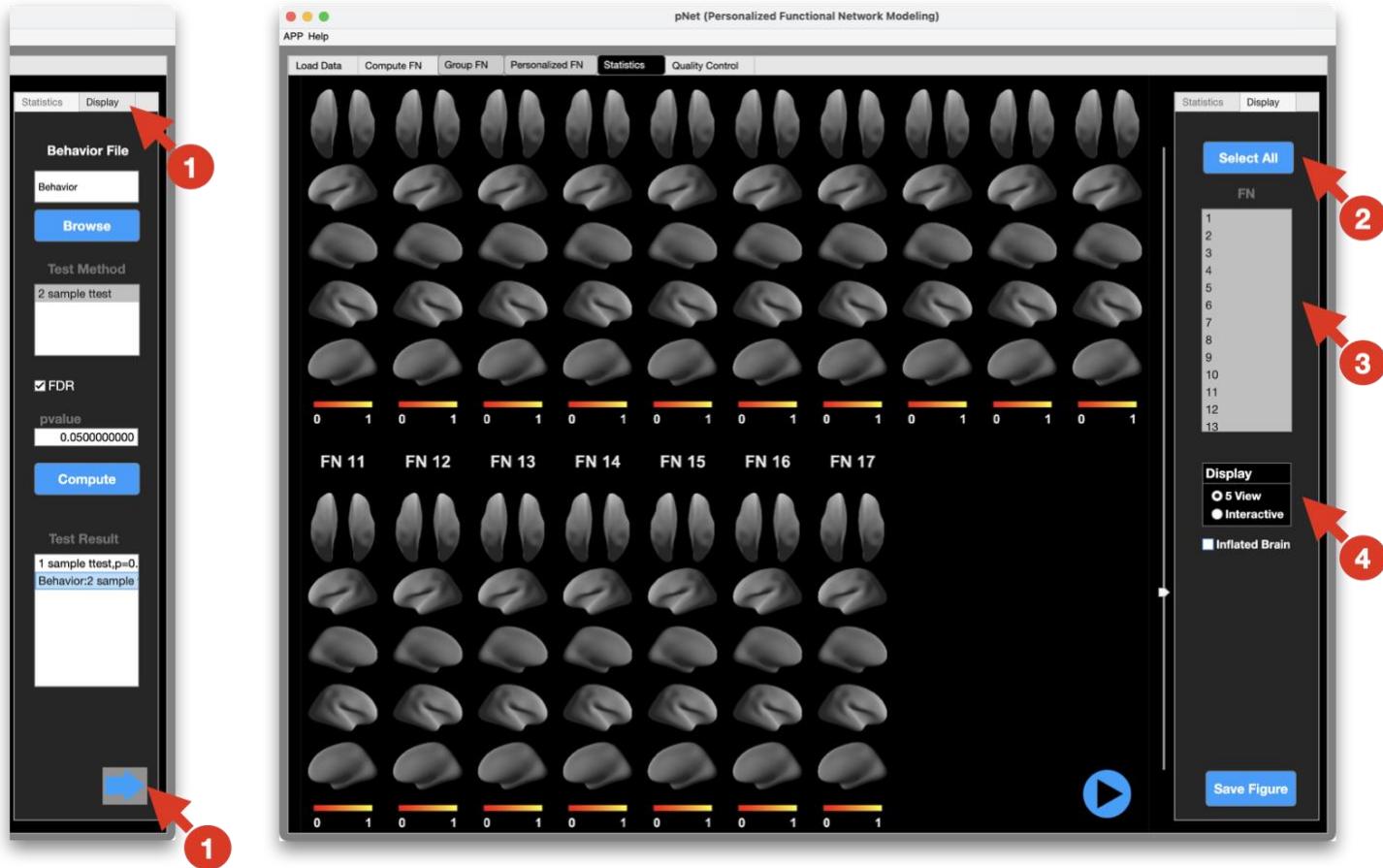


Check the two-sample t-test results

- (1) Navigate to the display module.
- (2) Click the select all button, if needed.
- (3) Select one or multiple FNs, if preferred.
- (4) Choose the display option.

Quality Control

The quality control is to ensure the spatial correspondence between personalized FNs and corresponding group-level FNs. Specifically, the personalized FN needs to have the highest spatial similarity (defined by Pearson correlation) to its corresponding group-level FN. The minimum spatial similarity value for all personalized FNs is noted as the $\text{min}\Delta\text{Sim}$ for each subject. If $\text{min}\Delta\text{Sim}$ is smaller or equal to 0, it means that there is at least one personalized FN mismatches to the group-level FN or has no higher spatial similarity.

Since our personalized FN computation method has quality assurance built in, all results will meet the quality control. Here we show one example with two mismatched FNs which will not be obtained from the toolbox.

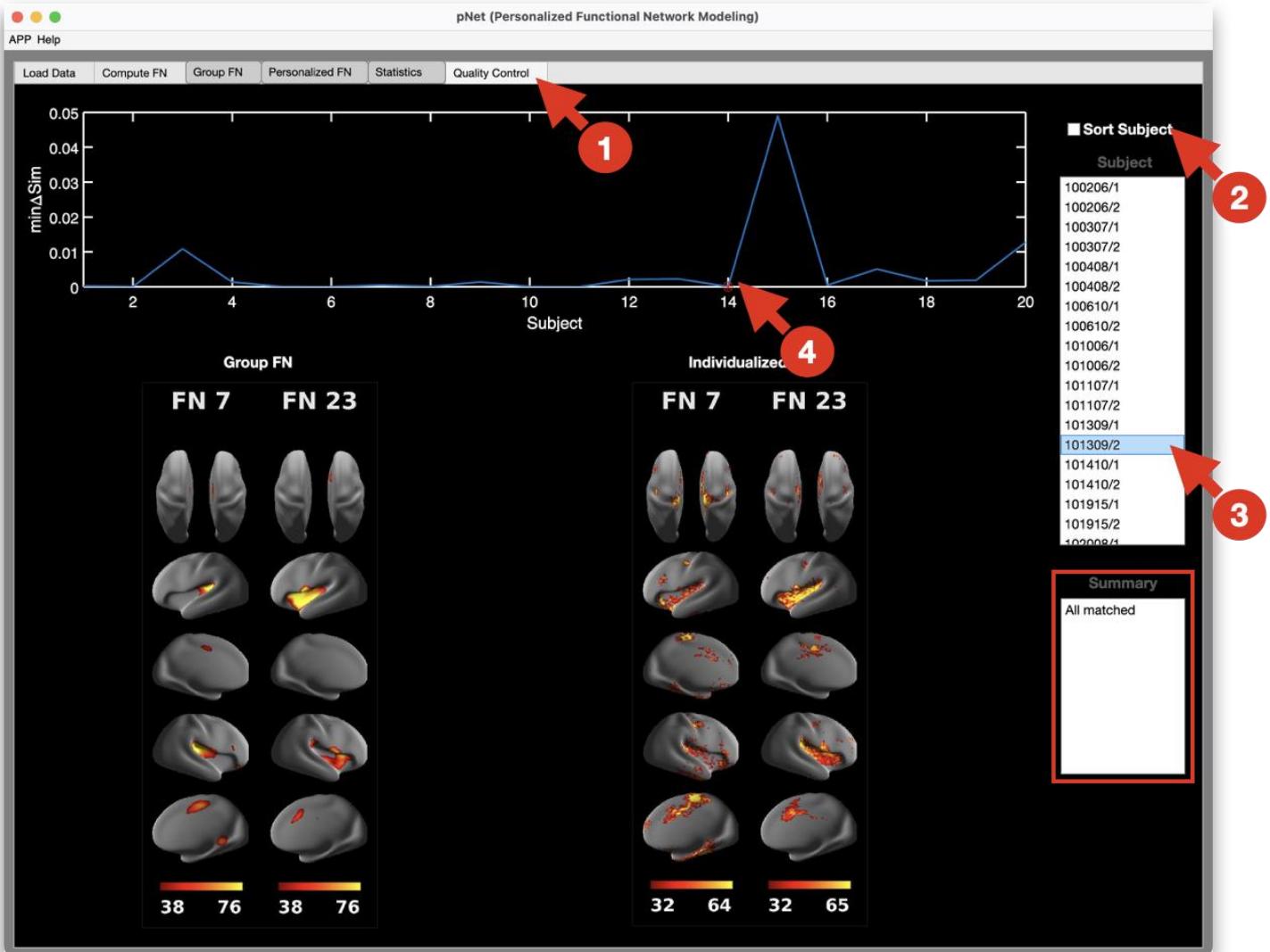

Check quality control.


- (1) Navigate to the quality control module.
- (2) Unselect the “Sort Subject” box, if needed. This allows to show quality control results of all subjects.
- (3) Choose subject to see the quality control result.
- (4) Alternatively, click on the top middle figure, quality control index ($\text{min}\Delta\text{Sim}$), to see navigate to the result of that subject.

Quality control summary of the selected subject is shown in the red box.

Example

We provided three examples for users to learn how to run the personalized FN computation, and navigate through the precomputed results. For each example, we provides examples of simulated fMRI data, and precomputed results of real fMRI data. Examples can be accessed on Google Drive: https://drive.google.com/drive/folders/1xkCy-0WqYvPA9ooq8txTdc0GsGC-YXMq?usp=share_link.

The data organization is as below.

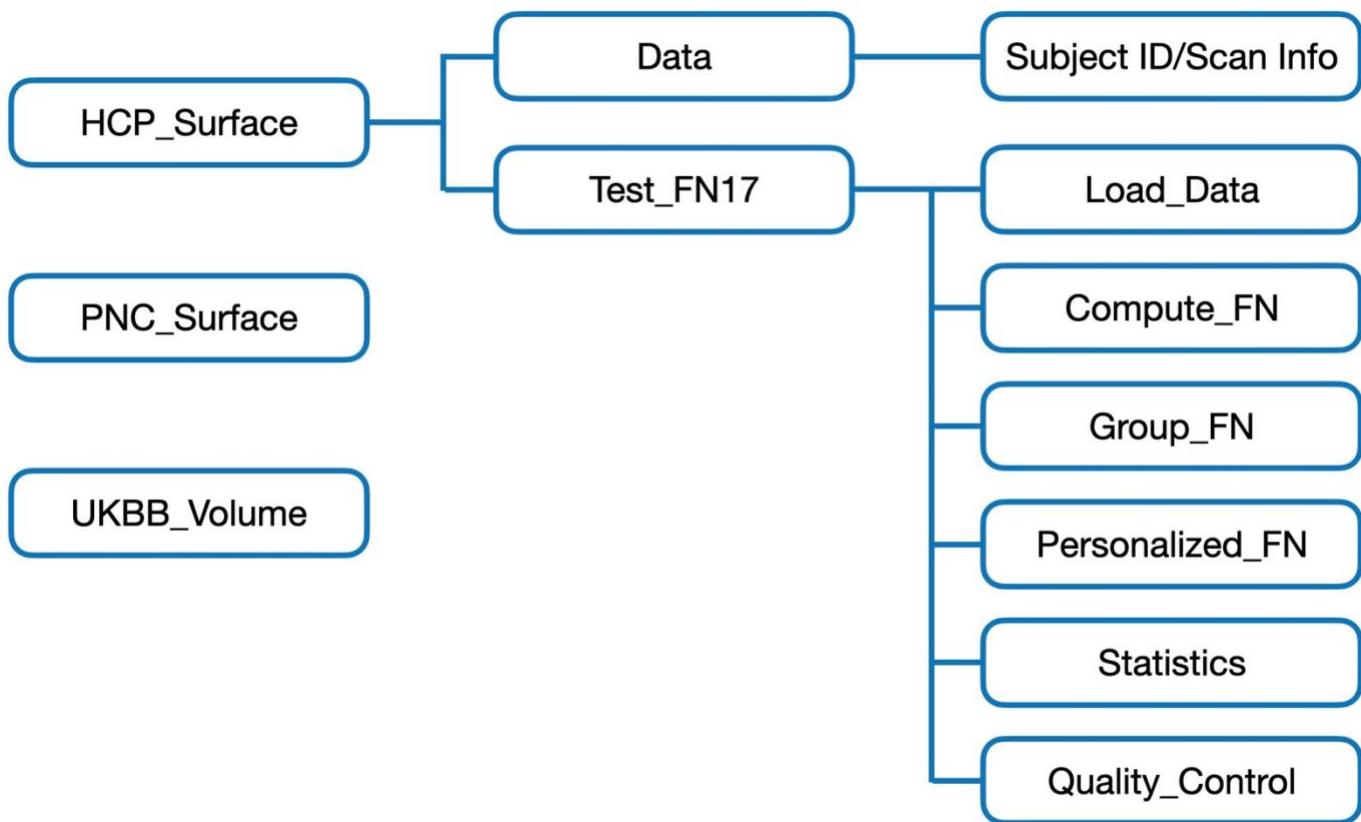

Data organization of examples

1. HCP surface data. (2 subject, 2 scans per subject, 400 volumes per scan) in .mat format.
Precomputed results of 10 subject data are stored in subfolder "Test_FN17".
2. PNC surface data (2 subject, 1 scan per subject, 300 volumes per scan) in .mgh format.
Precomputed results of 10 subject data are stored in subfolder "Test_FN17".
3. UKBB volume data (2 subject, 1 scan per subject, 200 volumes per scan) in .nii.gz format.
Precomputed results of 10 subject data are stored in subfolder "Test_FN17".

Each data folder contains simulated data organized by their subject ID, noted by a number, and scan information for multiple fMRI sessions or scans. For example,

Each precomputed result folder contains six subfolders. Each subfolder corresponds to one module in the pNet APP. For example, subfolder “Load_Data” contains setting information, fMRI scan list, brain mask file and others, which is required to be set in the data loading module.